



# ANTHEM<sup>®</sup>

## Distal Radius Fracture System



*Our mission is to deliver cutting-edge technology, research, and innovative solutions to promote healing in patients with musculoskeletal disorders.*

***Life moves us*** 

The Surgical Technique shown is for illustrative purposes only. The technique(s) actually employed in each case always depends on the medical judgment of the surgeon exercised before and during surgery as to the best mode of treatment for each patient. Additionally, as instruments may occasionally be updated, the instruments depicted in this Surgical Technique may not be exactly the same as the instruments currently available. Please consult with your sales representative or contact Globus directly for more information.

# SURGICAL TECHNIQUE GUIDE

## ANTHEM<sup>®</sup>

### Distal Radius Fracture System

System Overview .....	4
Implant Overview .....	9
Screw Compatibility .....	10
Volar Plating: Surgical Technique	
1. Approach .....	14
2. Retraction .....	16
3. Dissection .....	17
4. Fracture Reduction .....	20
5. Plate Application .....	20
6. Distal Provisional Fixation .....	22
7. Distal Fixation .....	24
Drill Guide Options .....	25
8. Proximal Fixation .....	26
9. Confirm Reconstruction .....	27
Final Construct .....	28
Optional: Removal .....	28
Lunate Facet Hook Plating: Surgical Technique	
Option A: Independent Facet Hook .....	31
Option B: Prior to Volar Plate Insertion .....	32
Option C: Following Volar Plate Insertion .....	33
Dorsal Plating: Surgical Technique .....	35
Lateral Plating: Surgical Technique .....	39
Bridge Plating: Surgical Technique .....	44
Ulnar Plating: Surgical Technique .....	49
Instrument Overview .....	54
ANTHEM <sup>®</sup> SS Distal Radius Implants .....	60
ANTHEM <sup>®</sup> Ti Distal Radius Implants .....	64
ANTHEM <sup>®</sup> Distal Radius Implants .....	66
ANTHEM <sup>®</sup> Distal Radius Fracture System Instruments .....	70
Important Information .....	74

# ANTHEM<sup>®</sup>

## Distal Radius Fracture System

The ANTHEM<sup>®</sup> Distal Radius Fracture System is a comprehensive fixation system designed to treat a variety of wrist fractures with anatomically contoured plates for intraoperative versatility.





## Plate Design

- Polyaxial holes offer a  $\pm 20^\circ$  cone of angulation
- ANTHEM® 7 Volar Plate designed for visibility of the intermediate column exposure
- ANTHEM® Double Row Volar Plate features an additional screw hole for maximum buttressing



## Unique Instruments

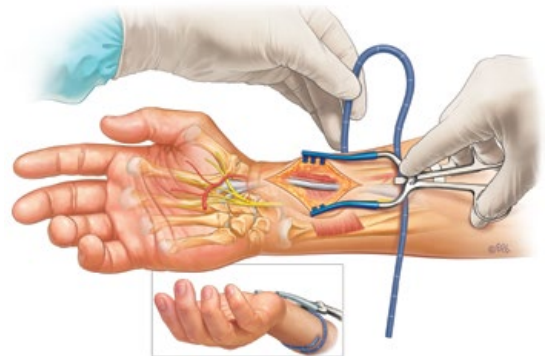
- One drill and driver to streamline procedure
- Stabilizing Radiolucent Weitlaner designed for secure retraction during the procedure
- Speed Lock Drill Guide quickly and rigidly attaches to the plate without threading



**T8 Self-Retaining Driver and 1.8mm Drill Bit**




**Speed Lock Drill Guide**



**Stabilizing Radiolucent Weitlaner**

## Optimal Case Flow

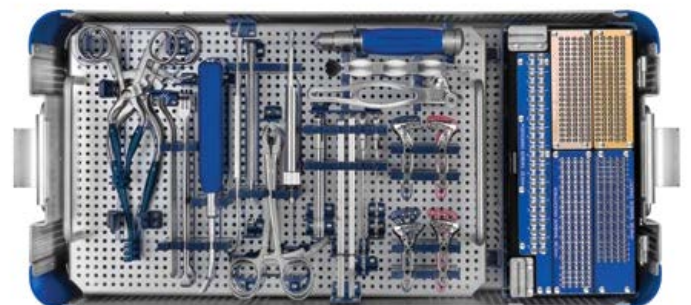- Graphic case trays are designed to optimize procedure flow



# Double Row Volar Plate

Anatomic plate features an additional screw hole to maximize buttressing.



## Maximum Buttressing



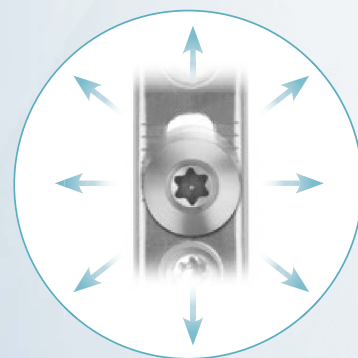
## Polyaxial Design

$\pm 20^\circ$  cone of angulation  
Allows for screw removal up to four times




## Positioning Slot and Screw

Allows multidirectional fine-tuning of plate position




## K-wire Hole Trajectories

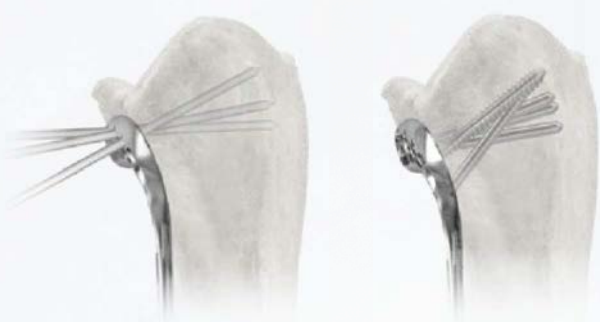K-wire holes match screw trajectories



## Radial Styloid

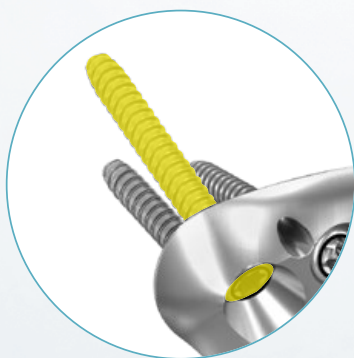Screw trajectory optimized for maximum purchase



# 7 Volar Plate

Anatomic single row volar plate designed to help increase intermediate column exposure of the distal radius.



## Innovative Plate Design

Allows clear view of fracture line



# Comprehensive Plate Portfolio



**Lunate Facet Hook Plate**



**Dorsal Plates**



**Lateral Plate**



**Bridge Plate**



**Ulna Plate**

# IMPLANT OVERVIEW

## Volar Plates

### **ANTHEM® 7 Volar Plate**

- Designed to help increase intermediate column exposure of the distal radius, allowing for a clear view of the fracture line.

### **ANTHEM® Double Row Volar Plate**

- Designed for maximized buttressing and includes a graft hole.



## Fragment Specific Plates

- Lunate Facet Hook Plate
- Dorsal (oblique and acute)
- Lateral
- Bridge
- Ulna



## Screws

- 1.8mm Locking Pegs
- 2.0mm and 2.5mm Locking Screws
- 2.5mm Non-Locking Screws
- 2.5mm Positioning Screws



Locking  
Peg

Locking  
Screws

Non-Locking  
Screw

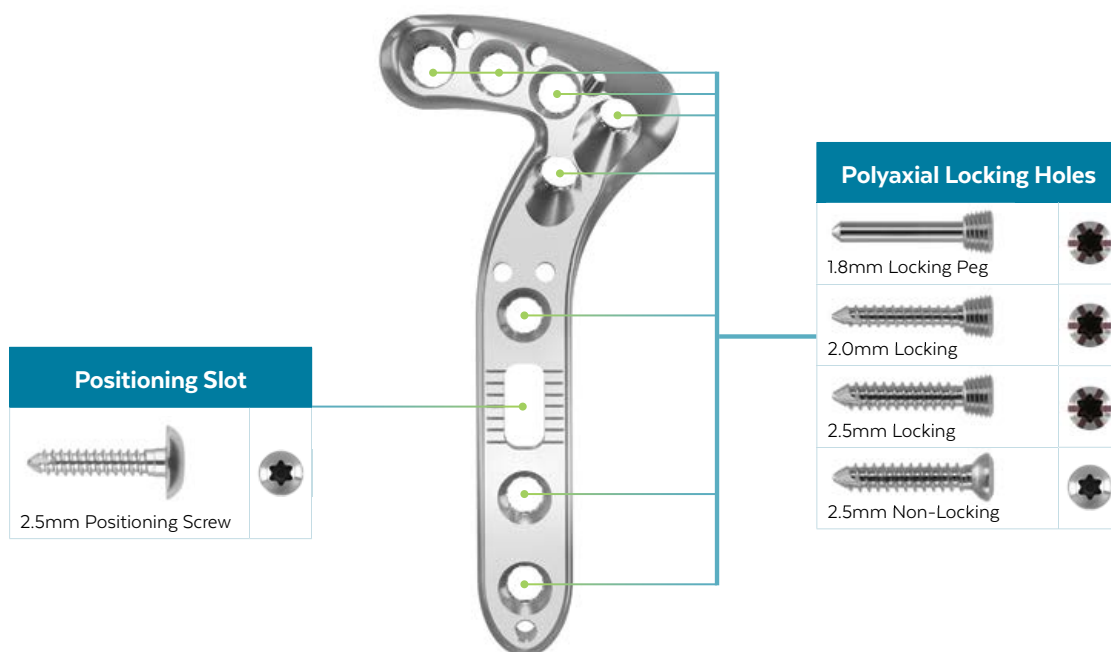
Positioning  
Screw




## SCREW COMPATIBILITY

If screw/plate locking is desired in a Polyaxial Locking Hole, only use Locking screws. MonoAx™ Locking Screws may not be used in a Polyaxial Locking Hole.

### ANTHEM® 7 Volar Plate




### ANTHEM® Double Row Volar

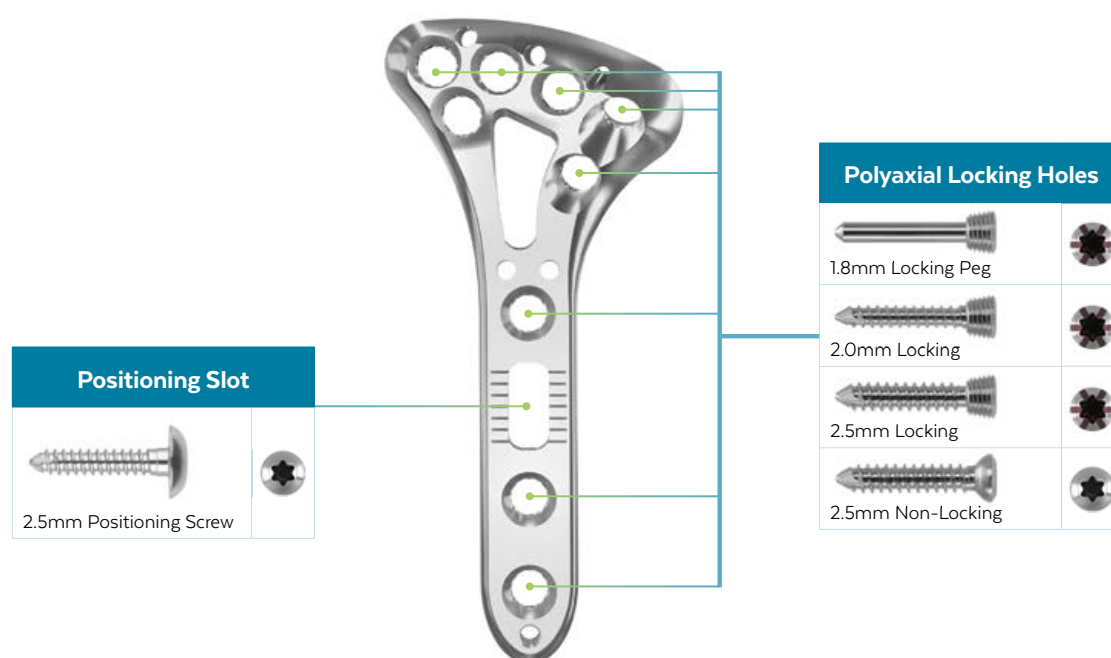


# SCREW COMPATIBILITY

If screw/plate locking is desired in a Polyaxial Locking Hole, only use Locking screws. MonoAx™ Locking Screws may not be used in a Polyaxial Locking Hole.
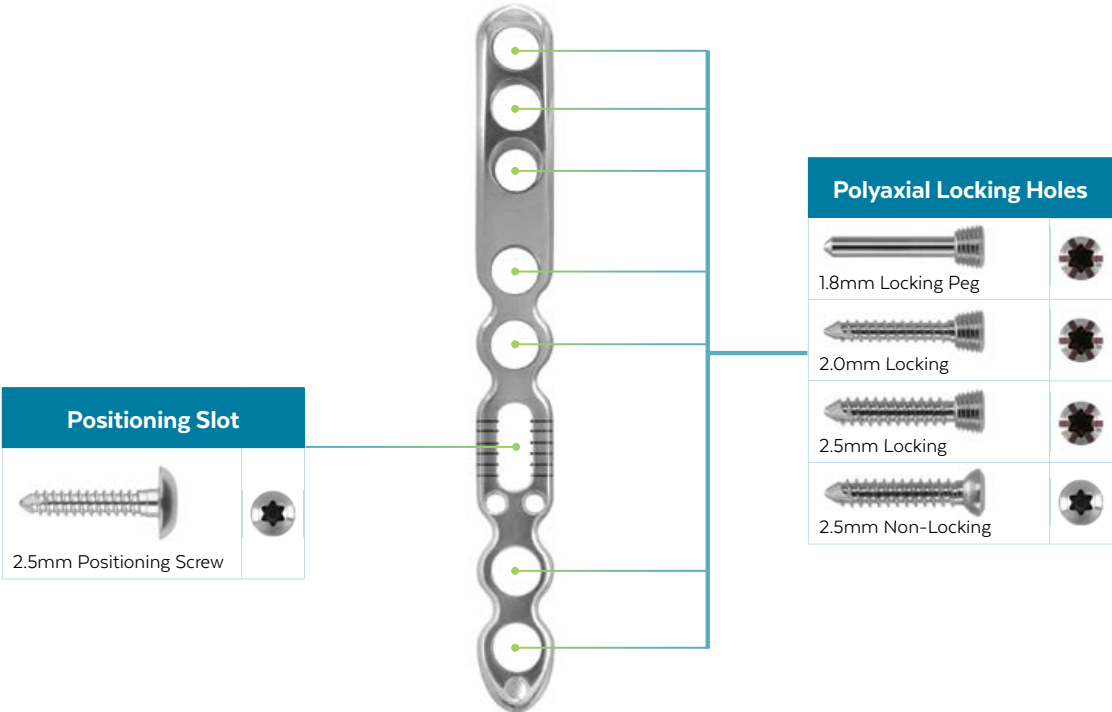
## Lateral Plate



## Dorsal Plates

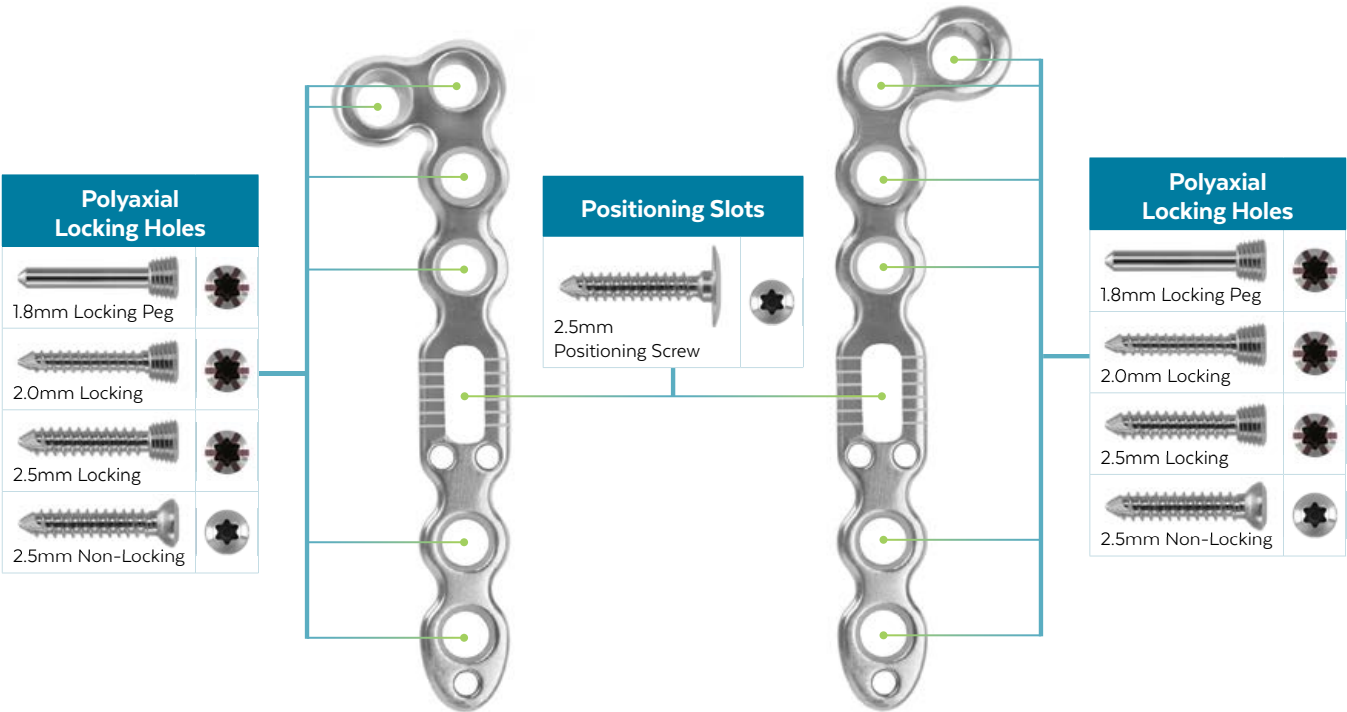


# SCREW COMPATIBILITY

If screw/plate locking is desired in a Polyaxial Locking Hole, only use Locking screws. MonoAx™ Locking Screws may not be used in a Polyaxial Locking Hole.

## Lunate Facet Hook Plate


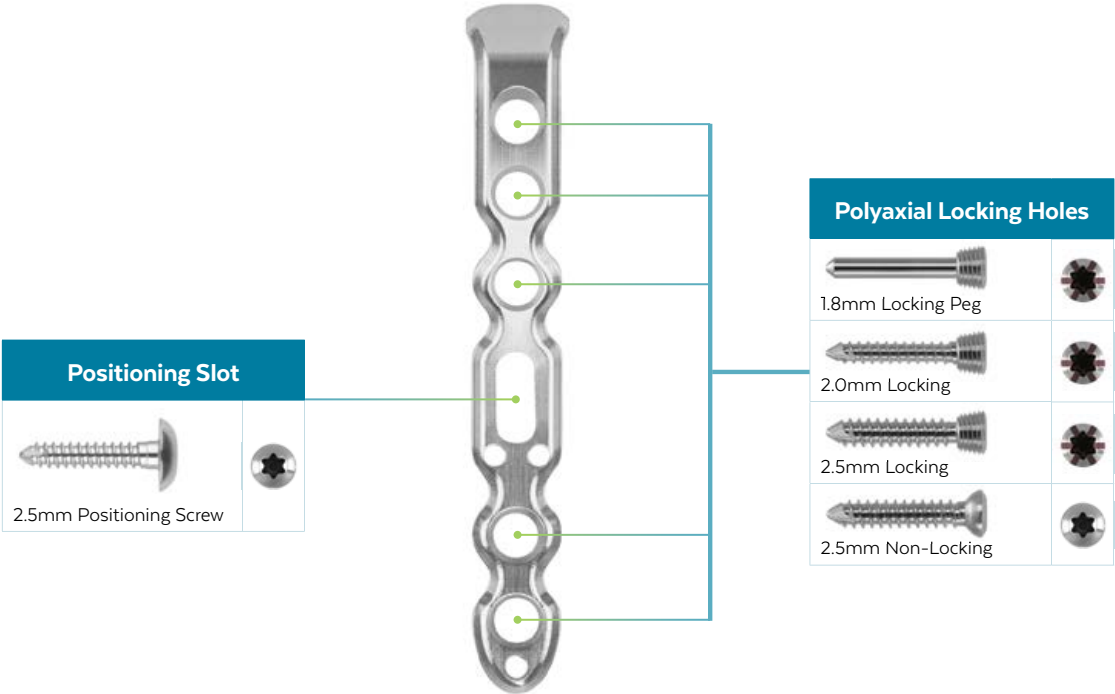
## Ulna Plate







# SCREW COMPATIBILITY

All holes in the Bridge Plate are MonoAx™ Locking Holes. If locking is desired, a MonoAx™ Locking Screw and/or a locking screw can be used.

## Bridge Plate



MonoAx™ Locking Holes	
 2.0mm Locking	
 2.5mm Locking	
 2.5mm Non-Locking	
 2.5mm MonoAx™ Locking Screw	

### OPTIONAL: SCREW INSERTION WITH 1.2Nm TORQUE LIMITER

The 1.2Nm Torque Limiter may be used to insert locking screws under power or in dense bone to help ensure proper tightening torque is not exceeded. Attach the T8 Driver to the 1.2Nm Torque Limiter. Under power, insert the locking screw until the maximum torque has been reached and an audible click is heard. Perform final tightening manually.

## SURGICAL TECHNIQUE

# ANTHEM<sup>®</sup>

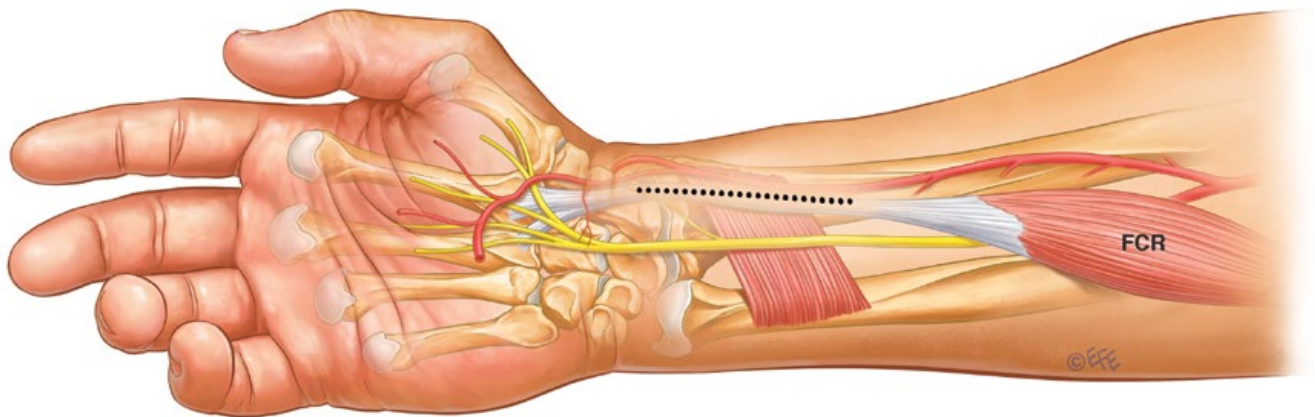
## Volar Plating

Refer to the package insert (also printed in the back of this manual) for important information on the intended use/indications, device description, contraindications, precautions, warnings, and potential risks associated with this system.

### STEP 1 APPROACH

The patient's arm is placed in a supine position.

Create an incision over the flexor carpi radialis (FCR) tendon using the FCR approach.



If necessary to cross the wrist, the incision should cross obliquely.



## ASSEMBLY OF THE STABILIZING RADIOLOUCENT WEITLANER

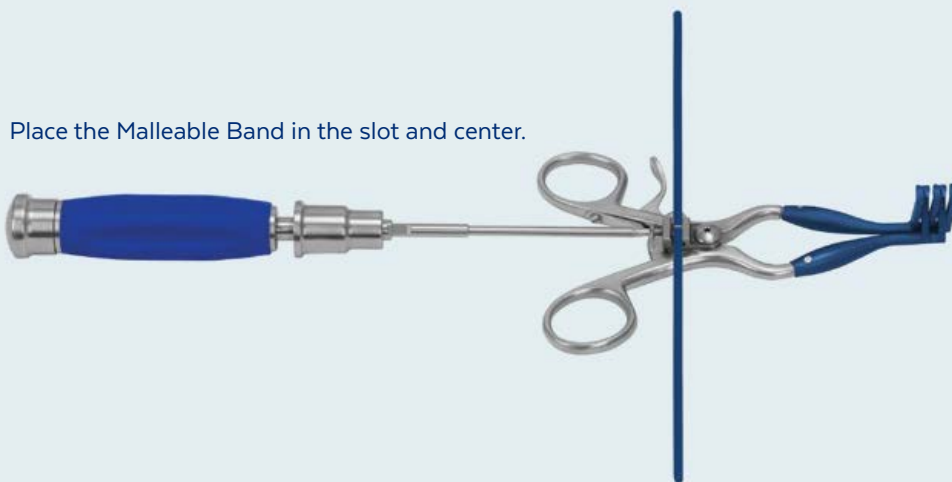
The Malleable Band provides intraoperative stability to the Weitlaner. Radiolucent arms eliminate the need for removal during fluoroscopy.

**Instruments:** T8 Driver and Malleable Band

Using the T8 Driver, loosen the set screw.



Place the Malleable Band in the slot and center.




Tighten the set screw with the T8 Driver.



Replacement bands are located on the middle tray.

## STEP 2 RETRACTION

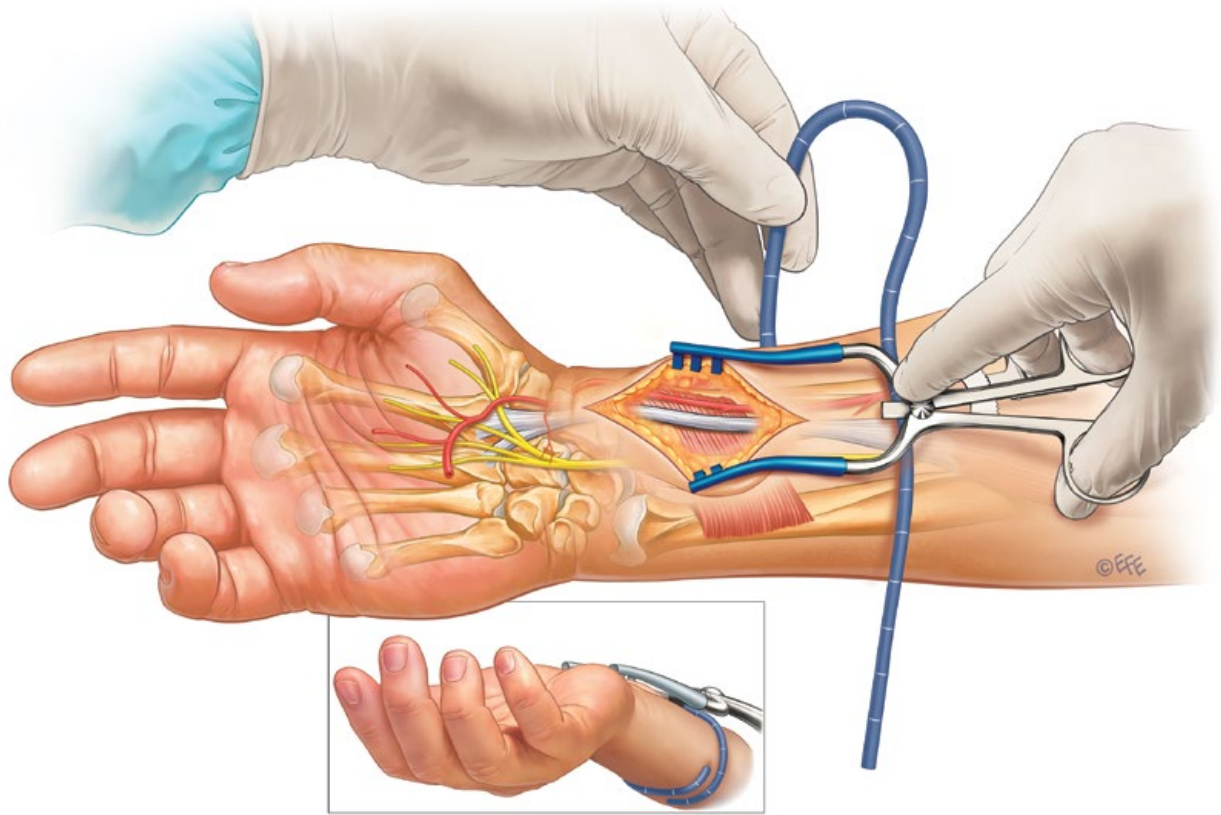Position the Stabilizing Radiolucent Weitlaner and retract. Wrap the Malleable Band around the patient's arm to secure the Weitlaner.




## STEP 3 DISSECTION

Three levels of dissection are used to access the radius: superficial, mid-level, and deep.

### Superficial Dissection

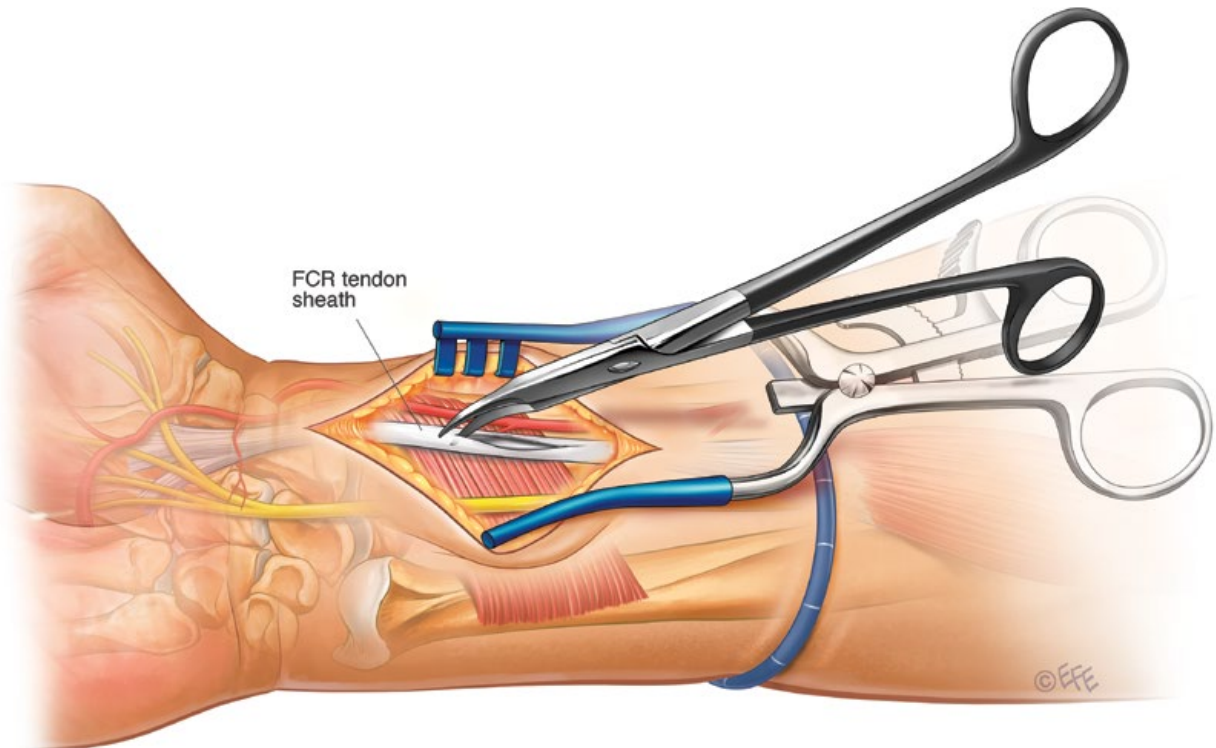Blunt dissection is used to access the FCR tendon. The FCR sheath is opened sharply.




Carefully avoid any aberrant branches of the palmar cutaneous branch of the median nerve crossing the FCR tendon sheath.

# DISSECTION (CONT'D)

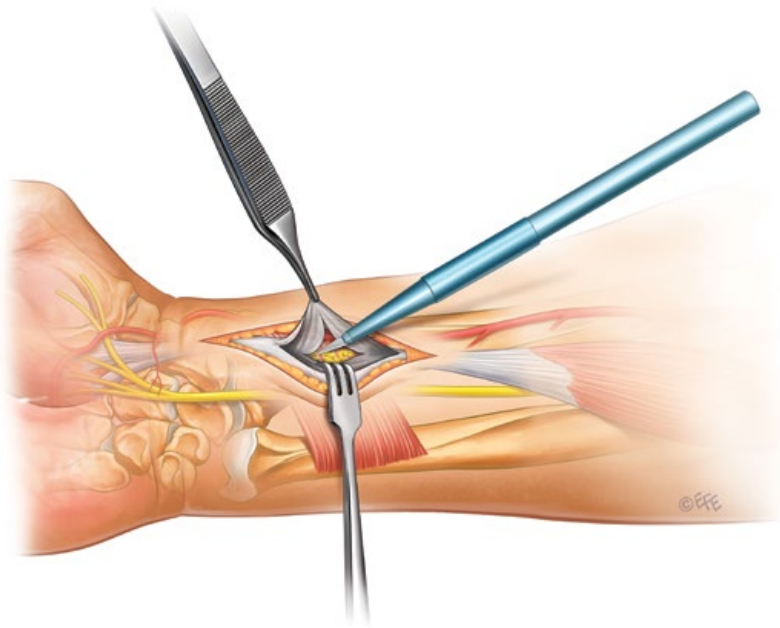
## Mid-Level Dissection

With the FCR tendon retracted, the floor of the tendon sheath is exposed and released sharply to reveal the deep volar compartment.



## Deep Dissection

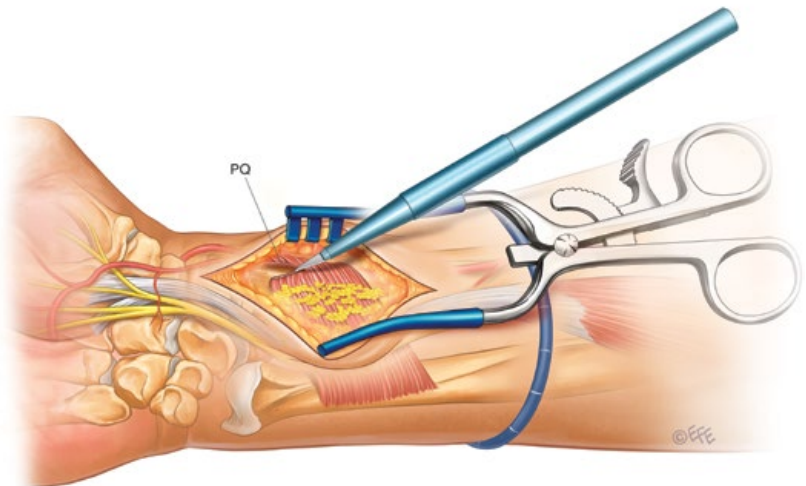
The plane between the flexor tendons and the pronator quadratus (PQ) muscle is bluntly developed. A loose fat pad may be encountered and excised. All flexor tendons, including the FCR tendon, should be retracted ulnarly and the radial artery retracted radially.



Exposure may be increased by releasing the FCR tendon sheath floor distally.

## Exposing the Fracture

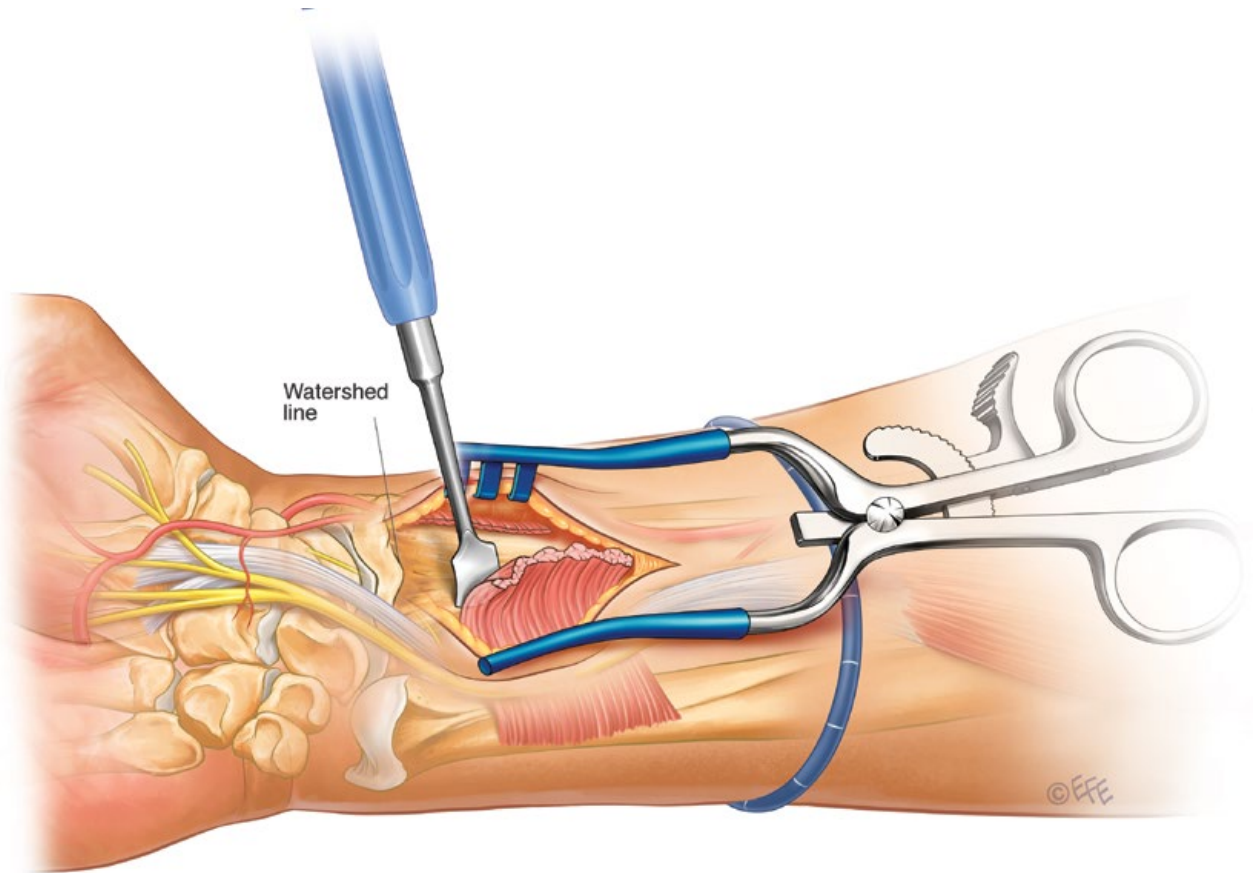The PQ muscle is released sharply across the radial border of the distal radius. Using the **Distal Radius Periosteal Elevator**, move the PQ muscle transversely, proximal to the watershed line.



Any remaining PQ muscle distal to the fracture line may be debrided away to improve visibility of the fracture line.



## STEP

## 4

## FRACTURE REDUCTION

Using fluoroscopy, reduce the fracture using the appropriate method for the fracture type. Ensure the radial length, alignment, and rotation are properly restored. Confirm reduction using fluoroscopy. K-wires may be used for temporary fixation.

## STEP

## 5

## PLATE APPLICATION

Select the plate that best accommodates the fracture pattern and patient anatomy. Using the **Plate Positioning Tool**, position the plate 2mm proximal to the watershed line on the radius.


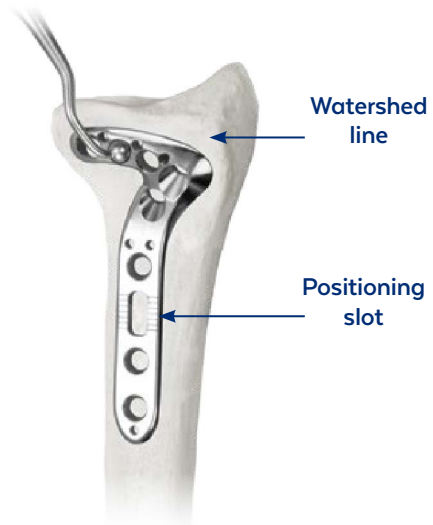


Plate position relative to watershed line

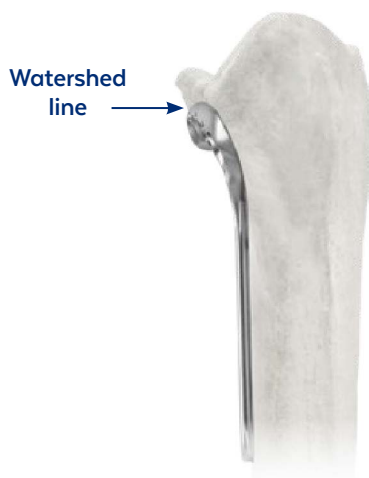Using the **1.8mm Drill Bit** and the **1.8mm Polyaxial Soft Tissue Protector**, drill a hole through the center of the positioning slot.



Drilling through fixed angle end of the 1.8mm Polyaxial Soft Tissue Protector into positioning slot

Use the **Depth Gauge** to determine the screw length.



Depth Gauge

Use the T8 Driver to select the desired 2.5mm Positioning Screw.

Use the **Measuring Gauge** to verify length.

Insert the positioning screw into the drilled hole using the T8 Driver. Verify plate placement using fluoroscopy. Adjust the plate proximal-distal or medial-lateral to the desired position. Tighten the positioning screw using the T8 Driver and 1.2Nm Torque-Limiting Handle. Once the plate is in an acceptable position with the positioning screw, the plate may be fixed definitively using additional shaft screws or provisionally held with K-wires through the plate shaft.

## POSITIONING SLOT AND SCREW

The positioning screw allows for intraoperative, multidirectional fine-tuning of the plate position.

The positioning hole and positioning screw allows for approximately 1.0mm of medial-lateral translation and 5mm of proximal-distal translation to optimize plate position.



Medial-lateral translation

## STEP

## 6

## DISTAL PROVISIONAL FIXATION

As an alternative to shaft-first fixation, the plate may be first applied to the distal fragment and then reduced to the radial shaft.

With the fracture reduced, the distal fragment is provisionally held with K-wires.

Using the K-wire attachment for the drill, insert the **1.6mm K-wire**. Insert the K-wire closest to the radial styloid. If desired, insert additional K-wires.



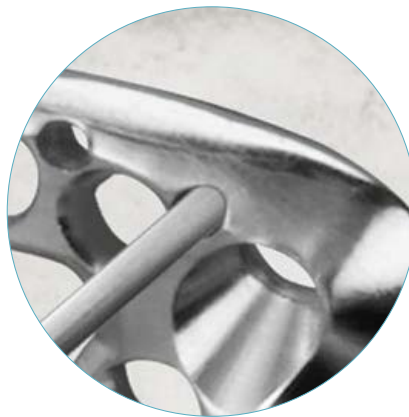
K-wire insertion

Confirm the trajectories of the K-wires using fluoroscopy.

The K-wire holes are fitted to match the nominal angle of subchondral screws. K-wires allow for preliminary fixation of the distal fragment and extrapolation of distal subchondral screw position.



K-wire trajectories



K-wire hole



If the fracture is fully mobilized prior to shaft plate application, the distal fragment may be reduced to the plate and provisionally fixed with K-wires before distal subchondral screw placement.

## STEP

## 7

## DISTAL FIXATION

Once the fracture is reduced, insert the distal subchondral screws by pre-drilling a hole using the 1.8mm Drill Bit and selected drill guide. Use the Depth Gauge to measure the desired screw length.

Insert 2.5mm Cortical Locking Screws, 2.0mm Locking Screws, or 1.8mm Smooth Locking Pegs using the T8 Driver. Verify screw and/or peg placement using fluoroscopy.

For final tightening, the 1.2Nm Torque-Limiting Handle is required.

If using a 2.0mm Locking Screw for anything other than buttressing, the **1.5mm Drill Bit** must be used.



Speed Lock Drill Guide



Depth Gauge



Prior to placement of additional subchondral locking screws, consider placing a non-locking cortical subchondral screw 2mm longer than measured and using as a reduction screw to help reduce the distal fragment and lag to the plate. This reduction screw may be exchanged for a locking screw. All subchondral locking screws are designed to be placed 2mm short of the dorsal surface of the distal radius to avoid tendon injury.



## DRILL GUIDE OPTIONS

### 1.8mm Speed Lock Drill Guide

This drill guide may be used to drill nominal trajectories. The thumb lock quickly locks the drill guide to the plate at the nominal trajectory.



Speed Lock Drill Guide



Speed Lock Drill Guide  
locked in place

### 1.8mm Polyaxial Soft Tissue Protector

This instrument allows for a 40° cone ( $\pm 20^\circ$ ) of angulation on the polyaxial end and nominal trajectory on the nominal end.



The angle should not exceed 20° from the central axis of the screw hole. Improper screw locking may result.



Polyaxial end



Fixed angle end

### Targeting Guide

The **Targeting Guide** may be used as a fixed angle drill and screw guide. Attach the corresponding Targeting Guide to the distal head. Use the **Targeting Guide Drill Sleeve** to drill each hole and insert screws through the Targeting Guide.



Targeting Guide



Targeting Guide Drill Sleeve

## STEP 8 PROXIMAL FIXATION

Shaft screw insertion varies depending on the reduction technique and fracture pattern. Pre-drill a hole using the 1.8mm Drill Bit and selected drill guide into the shaft holes. Measure hole depth using the Depth Gauge.

Insert the corresponding 2.0mm or 2.5mm Cortical Locking Screw or 2.5mm Non-Locking Screw using the T8 Driver.

For final tightening, the 1.2Nm Torque-Limiting Handle is required.

Verify placement using fluoroscopy.



1.8mm Drill Bit and  
Speed Lock Drill Guide



Screw inserted with  
the T8 Driver



The 1.8mm Drill Bit may be used in the proximal shaft holes and the distal subchondral holes. With the exception of the positioning slot, all screw holes accept locking or non-locking screws.

**STEP****9****CONFIRM RECONSTRUCTION**

Confirm joint reconstruction, screw placement, and trajectories using fluoroscopy. Use multiple views to verify proper reconstruction and confirm the distal screws have not penetrated the articular surface.



**Multiple views to confirm reconstruction**

## FINAL CONSTRUCT



## OPTIONAL: REMOVAL

If removal is required, use the T8 Driver to unlock the locking screws from the plate but do not remove the screws from the plate. This prevents simultaneous rotation of the plate. Use the T8 Driver to remove all non-locking and locking screws from the plate. Once all screws are removed, the plate can be removed.



Straight drivers are recommended for removal.

## SURGICAL TECHNIQUES

# ANTHEM<sup>®</sup>

## Fragment Specific Plating

The ANTHEM<sup>®</sup> Distal Radius Fracture System is a comprehensive system that includes fragment specific plates. These plates include the Lunate Facet Hook, Dorsal Acute and Dorsal Oblique, Lateral, Bridge, and Ulna plates. The surgical technique for each plate is provided in this section. Refer to the Volar plating technique for details on drill guides and removal instructions.



**Lunate Facet Hook Plate**



**Dorsal Plates**



**Lateral Plate**



**Bridge Plate**



**Ulna Plate**

## SURGICAL TECHNIQUE

# ANTHEM<sup>®</sup>

## Lunate Facet Hook Plating

### LUNATE FACET REPAIR

The Lunate Facet Hook Plate is designed for fracture patterns that involve the volar ulnar corner of the distal radius.

The plate internally fixes and/or buttresses the lunate facet with tines that hook over the rim of the facet.

#### Lunate Facet Hook Plate Application Methods:

- A) Independent Facet Hook
- B) Prior to Volar Plate Insertion
- C) Following Volar Plate Insertion

Lunate  
Facet Hook



**Option A**  
Independent Facet Hook



**Option B**  
Prior to Volar Plate Insertion



**Option C**  
Following Volar Plate Insertion



## Option A: Independent Facet Hook

This fragment specific technique is used for isolated fractures of the lunate facet and may be used independently to capture the facet.

### STEP 1 PLATE APPLICATION

With the fracture reduced, apply the plate using the **Lunate Facet Inserter**. Using the **Tamp** and **Mallet**, tamp the tines in or over the facet.

### STEP 2 PROXIMAL FIXATION

Drill a hole in the most proximal end of the slot in the plate using the Polyaxial Soft Tissue Protector and the 1.8mm Drill Bit. Measure the desired screw length with the Depth Gauge and insert a 2.5mm Non-Locking Screw using the T8 Driver. Confirm screw placement using fluoroscopy.

For final tightening, the 1.2Nm Torque-Limiting Handle is required.



Non-locking screws lag the plate to the bone and help to maintain fracture compression, avoiding late translation of the hook plate distally.



Lunate Facet Hook Plate applied to lunate facet



Tines placed using Tamp and Mallet



Proximal fixation using the 1.8mm Drill Bit and fixed angle end of the 1.8mm Polyaxial Soft Tissue Protector



Screw placement



Final fluoroscopic image



Final construct

## Option B: Prior to Volar Plate Insertion

This method is used when lunate facet fixation is necessary prior to Volar Plate application.

### STEP 1 PLATE APPLICATION

Apply the plate using the Lunate Facet Inserter. Using the Tamp and Mallet, tamp the tines in or over the facet.

### STEP 2 PROXIMAL AND PRELIMINARY FIXATION

Insert the 1.6mm K-wire. Using the K-wire for provisional fixation, place the Volar Plate over the K-wire onto the Lunate Facet Hook Plate.



Volar Plate positioned over the K-wire

### STEP 3 VOLAR PLATE APPLICATION

Follow the Plate Application steps beginning on page 20.



Volar Plate positioned over Lunate Facet Hook Plate



When positioned correctly, the ulnar-most subchondral locking screw aligns with the oblong hole of the Lunate Facet Hook Plate.

## Option C: Following Volar Plate Insertion

When an unstable lunate facet fracture is identified following volar plate fixation, the Lunate Facet Hook Plate may be applied under the Volar Plate.

### STEP 1 PLATE APPLICATION

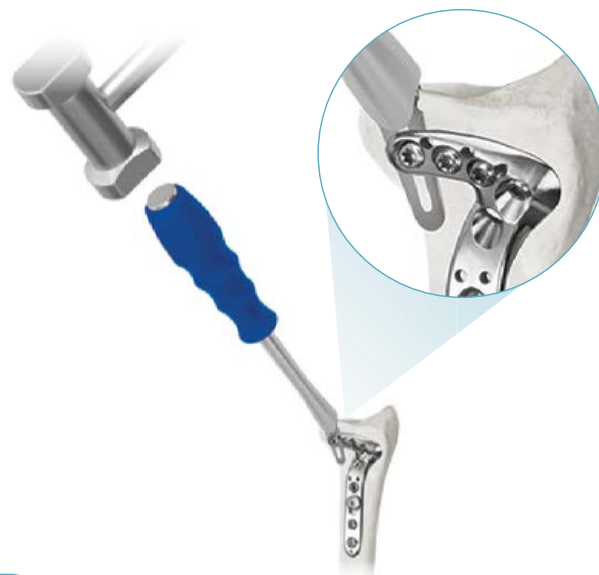
Apply the Lunate Facet Hook Plate under the Volar Plate, lining up the Lunate Facet Hook Plate with the ulnar-most subchondral screw hole.

Remove the ulnar-most subchondral Volar Plate locking screw using the T8 Driver.

Use the **Freer Elevator** to raise the distal end of the plate. Tamp the Lunate Facet Hook Plate under the Volar Plate using the Lunate Facet Inserter to direct the hook plate into position. Center the slot with the hole in the Volar Plate.



Ensure the Tamp does not contact the Volar Plate during impaction.



Tamping the Lunate Facet Hook Plate

### STEP 2 SCREW INSERTION

Use the T8 Driver to reinsert the ulnar-most subchondral 2.5mm Locking Screw into the Volar Plate.



Screw insertion

### STEP 3 CONFIRM PLACEMENT

Confirm placement of the Lunate Facet Hook Plate using fluoroscopy.



Final placement

## FINAL CONSTRUCTS



Independent Lunate Facet Hook Plate



Volar Plate and Lunate Facet Hook Plate

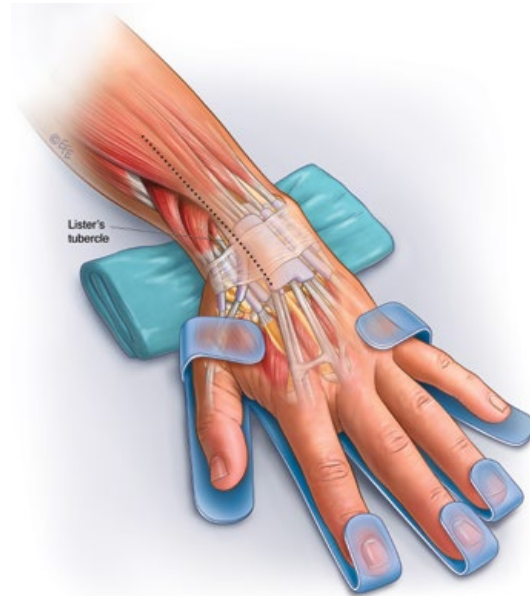
## SURGICAL TECHNIQUE

# ANTHEM<sup>®</sup>

## Dorsal Plating

### STEP 1 APPROACH

The patient's hand is positioned prone and secured with a lead hand or similar positioning device to maintain stability and orientation in the anterior-posterior (AP) plane.




Hand position

### STEP 2 EXPOSURE

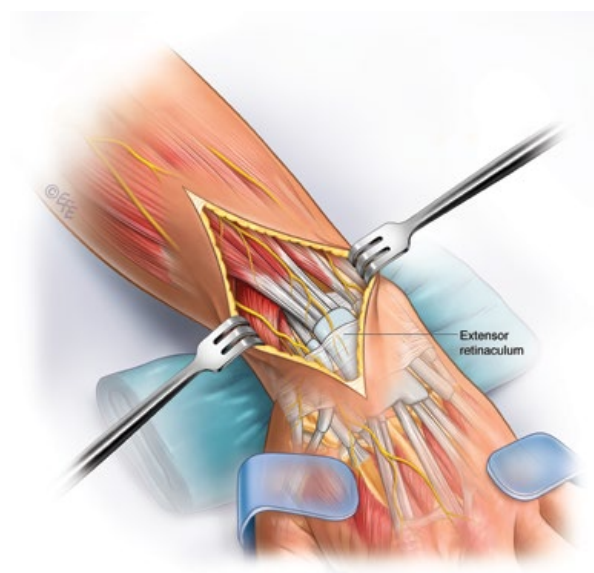Create a longitudinal incision approximately 8-10cm overlying Lister's tubercle.

Elevate full thickness skin flaps by carefully raising the soft tissue directly above the extensor retinaculum to preserve the dorsal radial and ulnar sensory nerves.




Retracting skin flaps  
to expose retinaculum

## EXPOSURE (CONT'D)

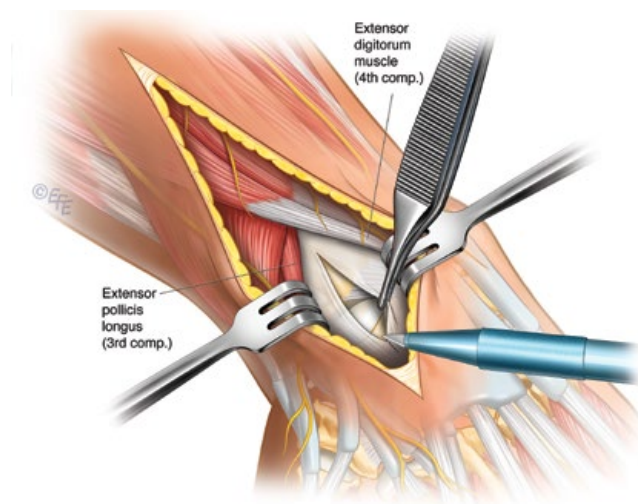Create a longitudinal incision between the third and fourth dorsal compartment or through the fourth compartment.



A double rectangular step cut allows for re-approximation of the retinaculum without transposition of the extensor pollicis longus (EPL) and allows tendon coverage.




Longitudinal incision

To view the articular surface, a dorsal capsulotomy may be performed.



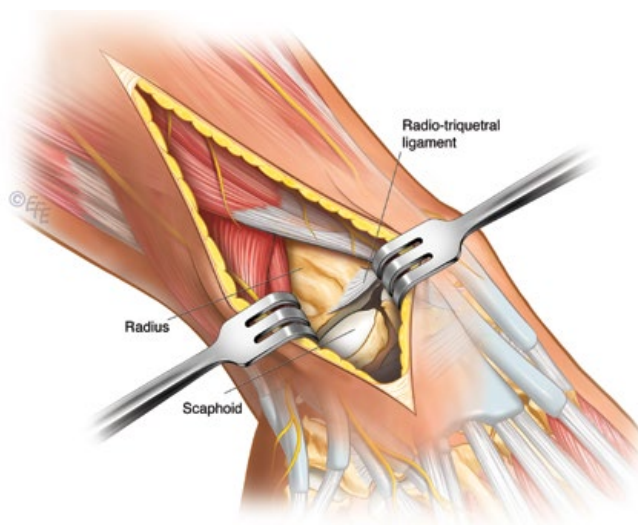Carefully preserve the dorsal radioulnar ligaments.



Wrist capsule

## STEP 3 REDUCTION

Using the **Dental Pick** or the Freer Elevator, elevate the articular fragments with allograft bone chips or suitable bone graft supplements.

Reduce the fragments under fluoroscopy using **Point-to-Point** or **Lobster Claw Reduction Forceps**. If preliminary fixation is required, use 1.6mm K-wires.



Traction on the wrist allows for direct visualization of the joint surface of the distal radius.



While reducing fragments under radiographic imaging, consider radial length (ulnar variance), articular reduction, palmar tilt, and radial inclination (height).



## STEP 4 PLATE SELECTION

Select an acute or oblique Dorsal Plate that best fits the patient's anatomy and fracture pattern.

## STEP 5 PLATE APPLICATION

Using the plate holding instrument, apply the plate to the midline of the dorsal and radial columns to buttress the fracture and maintain and restore volar and dorsal tilt.

If necessary, fix provisionally with K-wires and use fluoroscopy to determine correct placement.

Fix the plate through the positioning slot. Drill a hole in the center of the positioning slot using the 1.8mm Drill Bit and the 1.8mm Polyaxial Soft Tissue Protector. Using the Depth Gauge, measure hole depth. Insert the selected 2.5mm Positioning Screw with the T8 Driver. Adjust the plate proximal-distal or medial-lateral for optimal position. Tighten the positioning screw using the T8 Driver and the 1.2Nm Torque-Limiting Handle.

## STEP 6 PLATE FIXATION

### Proximal Fixation

Insert proximal screws by pre-drilling a hole using the 1.8mm Drill Bit and selected drill guide. Use the Depth Gauge to measure the desired screw length.

Insert a 2.5mm Cortical Non-Locking Screw or a 2.5mm Locking Screw using the T8 Driver.

### Distal Fixation

Insert distal subchondral screws by pre-drilling a hole using the 1.8mm Drill Bit and selected drill guide. Use the Depth Gauge to measure the desired screw length.

Insert 2.5mm Cortical Locking Screws, 2.0mm Locking Screws, or 1.8mm Smooth Locking Pegs using the T8 Driver. If using a 2.0mm Locking Screw for anything other than buttressing, the 1.5mm Drill Bit must be used.

Confirm placement using fluoroscopy and tighten all proximal and distal screws using the 1.2Nm Torque-Limiting Handle.



2.0mm Locking Screws may be used to capture smaller fragments distally.



Application of oblique Dorsal Plate



Pre-drilling of positioning slot using the Soft Tissue Protector



Depth Gauge measurement

## FINAL CONSTRUCT



Dorsal Plates

## SURGICAL TECHNIQUE

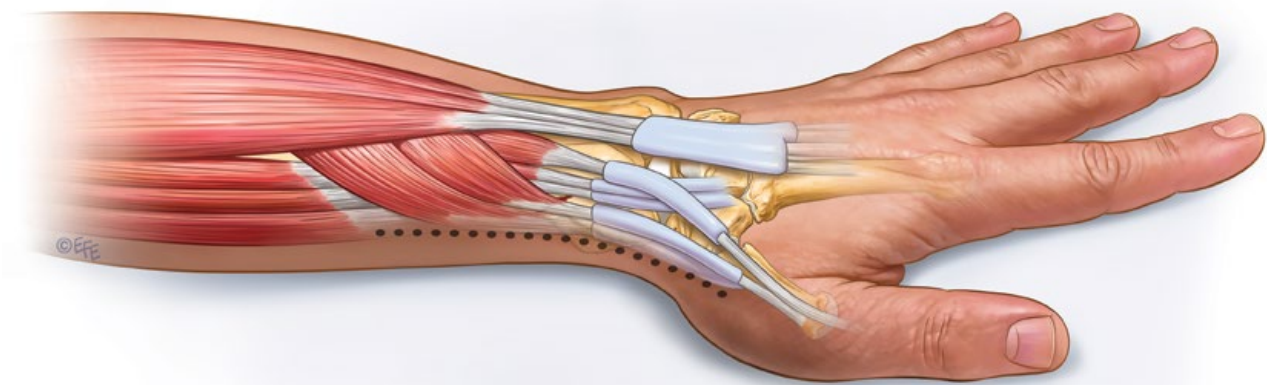
# ANTHEM<sup>®</sup>

## Lateral Plating

The Lateral Plate is designed for direct buttressing of distal radius fragments.

### STEP 1 APPROACH

Create a lateral incision directly along the radial border of the wrist.



If the Lateral Plate is applied with the Volar Plate, the Lateral Plate may be placed through the same volar incision.


## STEP 2 DISSECTION
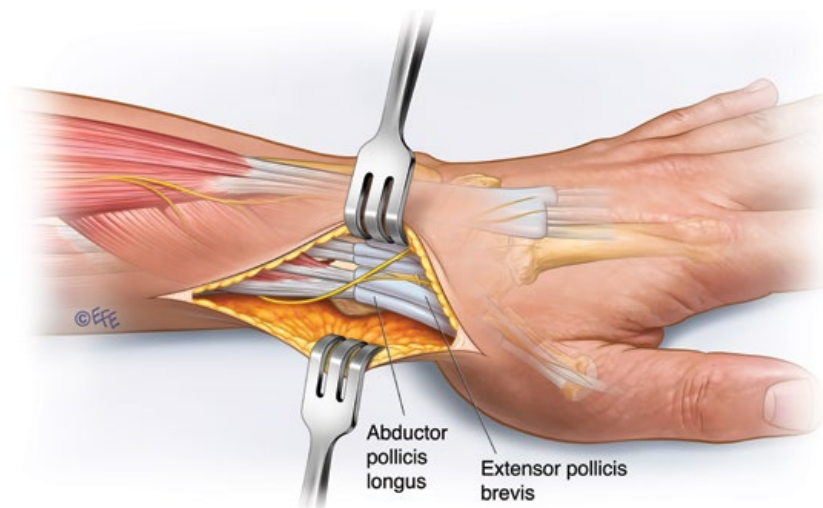
### Superficial Dissection

Blunt dissection is performed down to the first dorsal compartment of the wrist that contains the abductor pollicis longus (APL) and extensor pollicis brevis (EPB) tendons.




Branches of the radial sensory nerve are present in the surgical field and must be carefully dissected and protected.



### Deep Dissection

Release the first dorsal compartment and retract the APL and EPB tendons to expose the brachioradialis tendon.

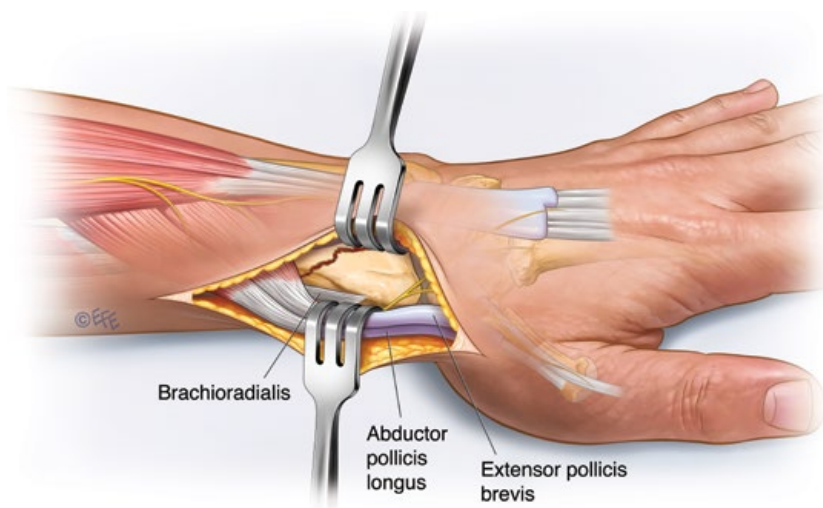The brachioradialis tendon may be released for fracture exposure, reduction, and plate placement.



## STEP 3 PLATE APPLICATION

Reduce the fracture. Apply the plate to the lateral border of the radial shaft.



Application of Lateral Plate

Using the 1.8mm Drill and selected drill guide, drill a hole through the positioning slot. Use the Depth Gauge to measure depth and select the desired 2.5mm Positioning Screw.



Pre-drilling of positioning slot

Insert the positioning screw using the T8 Driver. Adjust the plate proximal-distal or medial-lateral for optimal position.

The **Plate Positioning Tool** may be used to stabilize the plate during screw insertion. Confirm plate placement using fluoroscopy.



Positioning screw insertion

## STEP 4 PLATE FIXATION


### Proximal Fixation

Insert proximal screws by pre-drilling a hole using the 1.8mm Drill Bit and selected drill guide. Use the Depth Gauge to measure the desired screw length.

Insert 2.5mm Cortical Non-Locking Screws or 2.5mm Locking Screws using the T8 Driver.



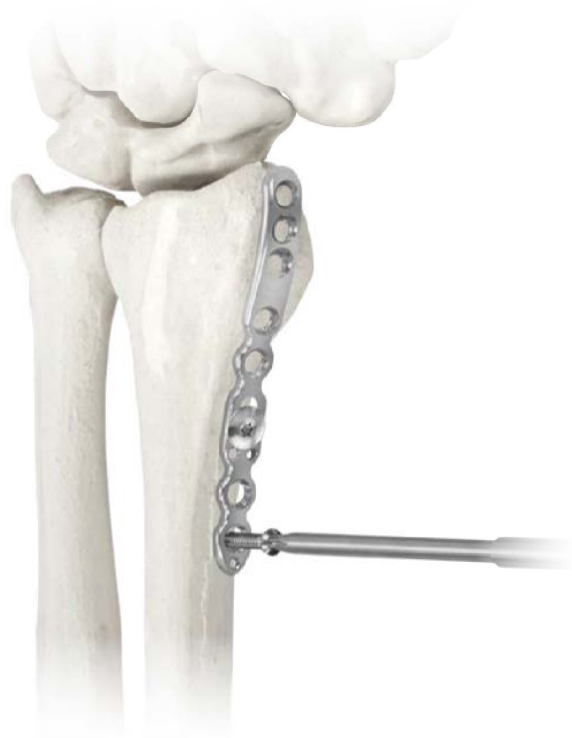Optimal plate position is achieved when the distal end of the plate meets the distal end of the radial styloid.



### Distal Fixation

Insert distal subchondral screws by pre-drilling a hole using the 1.8mm Drill Bit and selected drill guide. Use the Depth Gauge to measure the desired screw length.

Insert 2.5mm Cortical Locking Screws, 2.0mm Locking Screws, or 1.8mm Smooth Locking Pegs using the T8 Driver. If using a 2.0mm Locking Screw for anything other than buttressing, the 1.5mm Drill Bit must be used.

Confirm placement using fluoroscopy and tighten all proximal and distal screws.



The nominal angle for the distal subchondral screws is retrograde below the articular surface of the radiocarpal joint. For greater than nominal screw angulation, avoid inadvertent injury to the articular surface.





STEP

5

## CONFIRM REDUCTION

Using fluoroscopy, verify screw trajectories and confirm plate placement.



## FINAL CONSTRUCT



Lateral Plate

## SURGICAL TECHNIQUE

# ANTHEM<sup>®</sup>

## Bridge Plating

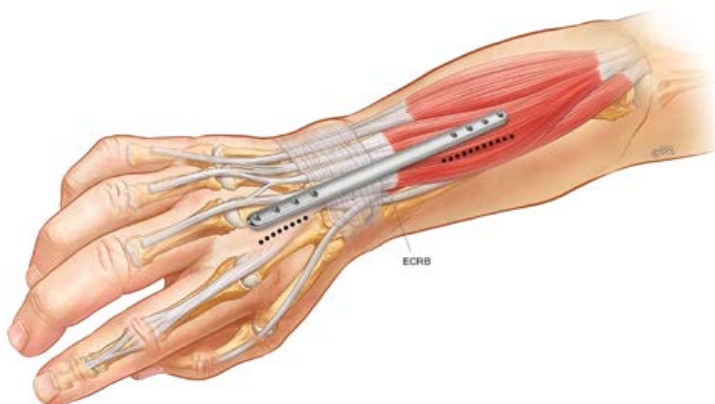
The Bridge Plate is designed for highly comminuted intra-articular fractures of the distal radius.

### STEP 1 APPROACH

Place the plate on the dorsal surface of the hand and forearm as a guide to placement. Mark the incision line using a surgical marker.



Use the index or middle metacarpal as a landmark. The incision should be made distally between the index and middle metacarpals to allow for access and placement.



Distal incision with Bridge Plate

### STEP 2 DISTAL EXPOSURE

Create a distal incision.

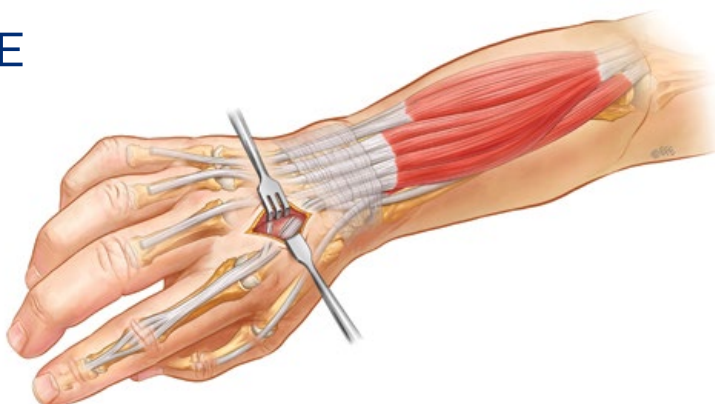
Retract the tendons to perform blunt dissection to the extensor tendons.



If the index metacarpal is selected for plate fixation, all extensor tendons may be retracted ulnarly.



If the middle metacarpal is selected, the extensor tendons to the index and middle fingers may require separation by dividing the juncturae tendinum. If released, the tendinum should be repaired upon closure to help avoid an extensor lag.




Distal incision and tendon retraction

## STEP


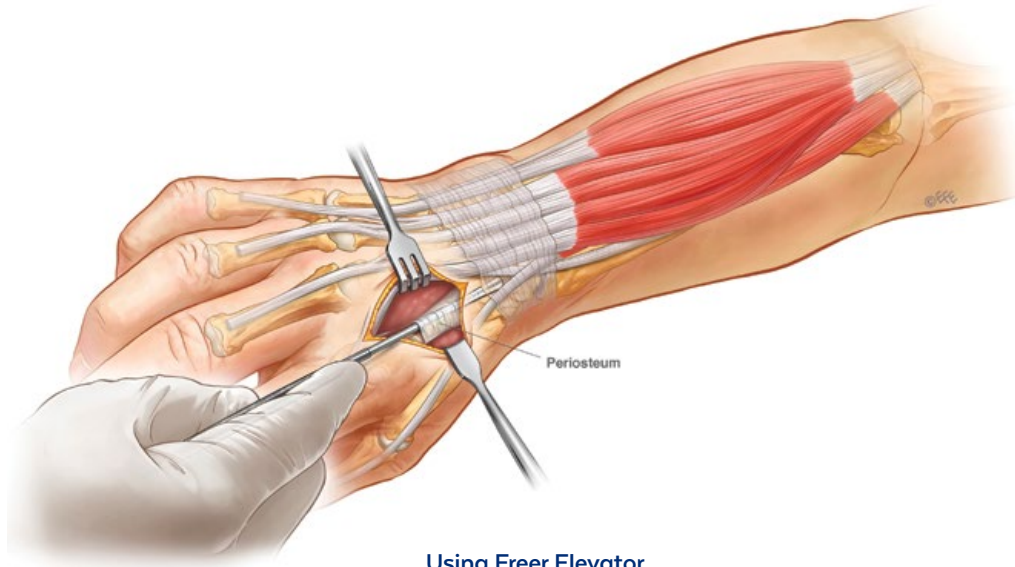
## 3

## PLATE APPLICATION

With the extensor tendons retracted, the metacarpal is exposed subperiosteally. Position the **Distal Radius Freer Elevator** retrograde under the extensor retinaculum to create a path for plate placement.



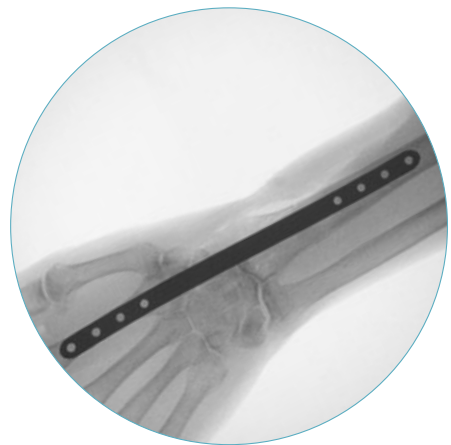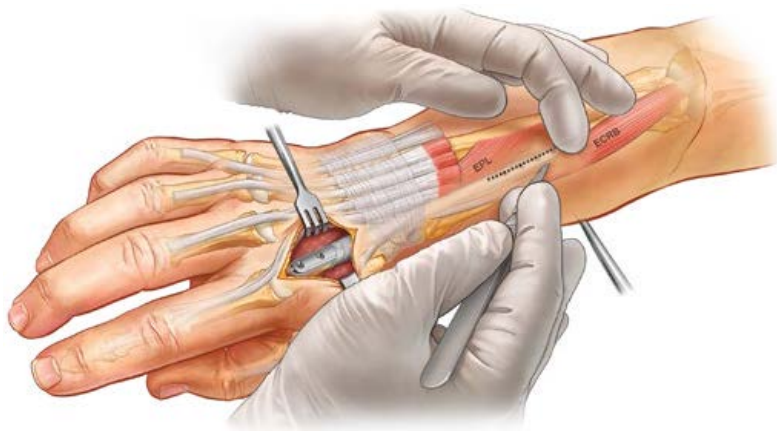Insert the Bridge Plate through the incision and align with the radius. Palpate the Bridge Plate on the radius for proximal placement. Confirm placement using fluoroscopy.



The Bridge Plate may be fixed to the index or middle metacarpal. Fixation to the index metacarpal allows for ease of exposure, tendon retraction, and access to the radial shaft through the second dorsal compartment. Wrist fixation in slight ulnar deviation may help to strengthen the power grip.

## STEP

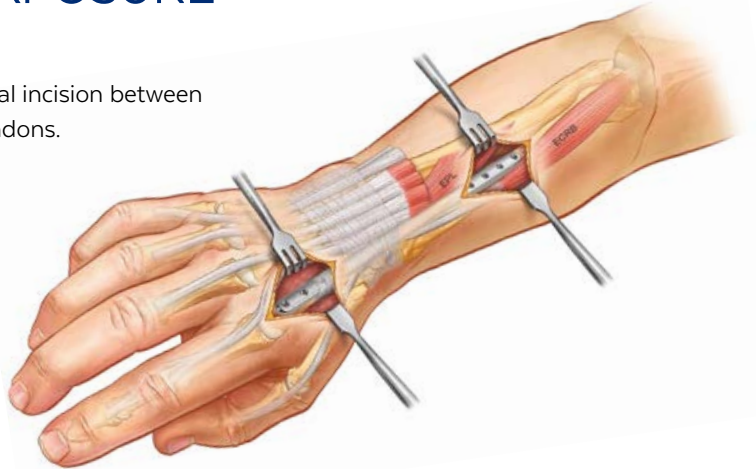
4

## PROXIMAL EXPOSURE

With retrograde placement satisfied, create a proximal incision between the extensor carpi radialis brevis (ECRB) and EPL tendons.



Place the Hohmann or Senn Retractor on either side of the radius to help center the plate.



Bridge Plate centered on radius

## STEP

5

## PLATE FIXATION

### Distal Fixation

Thread the 1.8mm Threaded Drill Guide into the second most distal monoaxial hole of the plate. Using the 1.8mm Drill Bit, drill through the second most distal hole of the plate. Measure hole depth using the Depth Gauge.

Using the T8 Driver, insert a 2.5mm Cortical Non-Locking Screw.

Confirm screw placement using fluoroscopy.



Distal hole fixation



Using the Depth Gauge



Prior to plate fixation, confirm satisfactory closed reduction of the fracture. If necessary, percutaneous or limited open reduction with K-wires or the Freer may be performed.



The holes in the Bridge Plate are monoaxial, which require the 1.8mm Threaded Drill Guide.

## Proximal Fixation

Thread the 1.8mm Threaded Drill Guide into the second most proximal monoaxial hole. Using the 1.8mm Drill Bit, drill through the second most proximal hole. Measure the desired screw length using the Depth Gauge.

Using the T8 Driver, insert a 2.5mm Cortical Non-Locking Screw.

Confirm screw placement using fluoroscopy.



Gentle traction and/or a reduction maneuver may be performed prior to proximal fixation.



Drilling second most proximal hole

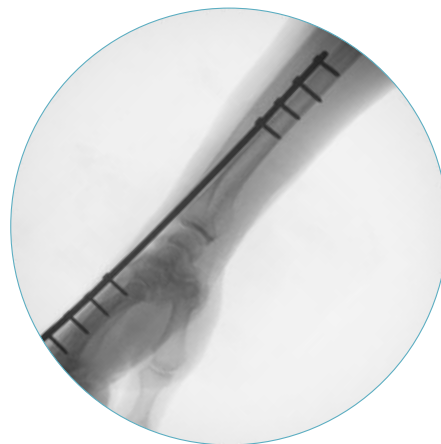
## STEP 6 SCREW INSERTION

Once fracture reduction and plate placement are achieved, fill additional screw holes proximally and distally with 2.5mm Locking Screws using the 1.8mm Drill Bit and the selected drill guide. Confirm screw placement using fluoroscopy.

For final tightening, the 1.2Nm Torque-Limiting Handle is required.



Locking screw fixation



Final placement



Following initial cortical screw fixation proximally and distally, fill the remaining holes with locking screws for stability.



# FINAL CONSTRUCT



Bridge Plate



## SURGICAL TECHNIQUE

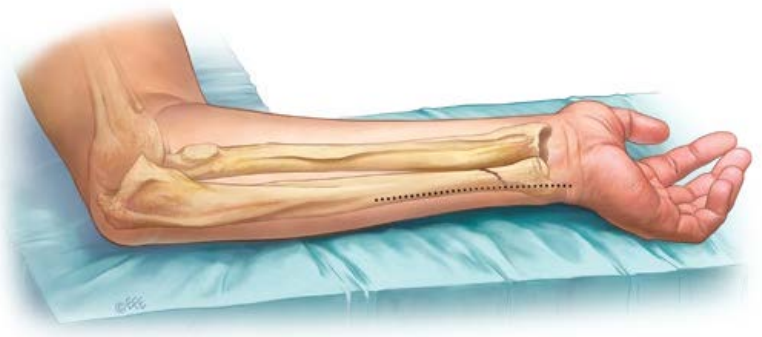
# ANTHEM<sup>®</sup>

## Ulnar Plating

Ulna Plates are used for fixation of an unstable ulna following distal radius repair.


### STEP 1 APPROACH

Using the subcutaneous ulnar approach, position the arm on a hand table and flex the elbow. Supinate the forearm to expose the subcutaneous border of the ulna.



### STEP 2 EXPOSURE

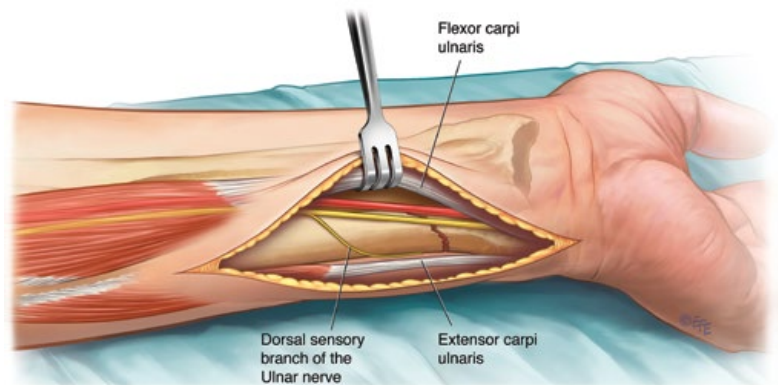Create a longitudinal incision approximately 1cm longer than the plate, distally and proximally. Split the interval between the extensor carpi ulnaris (ECU) and the flexor carpi ulnaris (FCU) to expose the ulnar shaft. The plate may be applied more dorsally, if desired.




Carefully avoid the dorsal sensory branch of the ulnar nerve.

**STEP****3****FRACTURE REDUCTION**

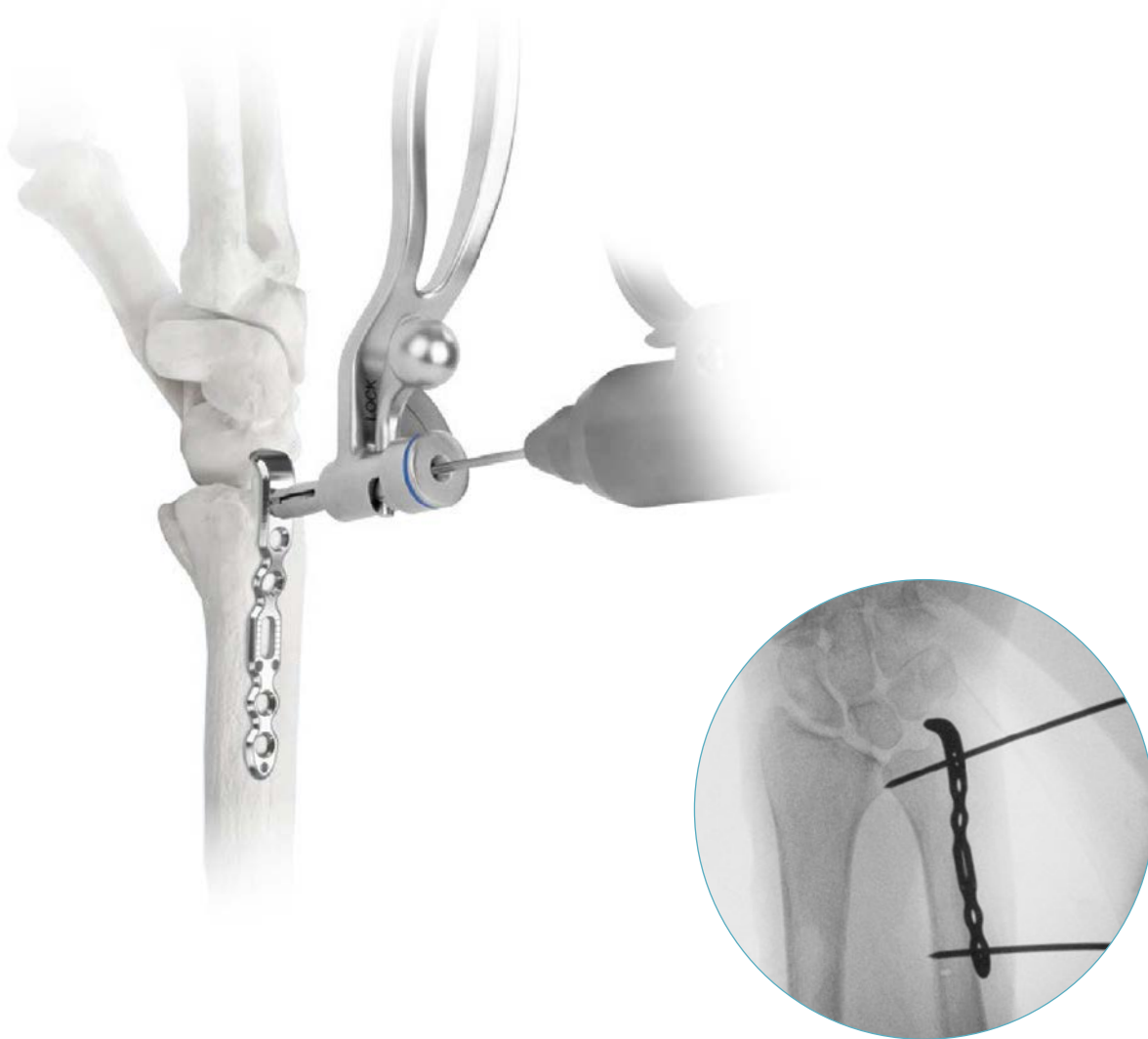Reduce the fracture and confirm reduction using fluoroscopy.

**STEP****4****PROVISIONAL PLATE FIXATION**

Using the Speed Lock Drill Guide, place 1.6mm K-wires in the most distal and proximal screw holes to provisionally hold the plate in position. Confirm K-wire placement using fluoroscopy.



**K-wire placement**



Using the most distal hole, a K-wire may be used for preliminary fixation.

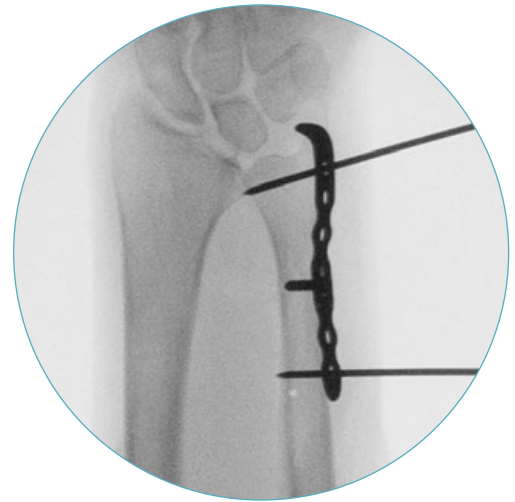
## STEP 5 PLATE APPLICATION

Using the 1.8mm Drill Bit and Soft Tissue Protector, drill a hole through the center of the positioning slot.

Using the Depth Gauge, measure the desired hole depth.

Insert the positioning screw using the T8 Driver. Adjust the plate proximal-distal or medial-lateral for optimal placement.

Confirm placement using fluoroscopy.




Positioning screw placement

## STEP 6 SCREW INSERTION

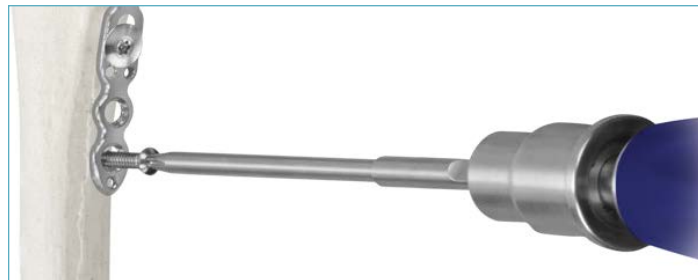Pre-drill for the remaining screws using the 1.8mm Drill Bit. Insert 2.0mm or 2.5mm Locking Screws. If using a 2.0mm Locking Screw for anything other than buttressing, the 1.5mm Drill Bit must be used.


For final tightening, the 1.2Nm Torque-Limiting Handle is required.

Confirm screw length and placement using fluoroscopy. Replace K-wires with locking screws.



**STEP****7****CONFIRM REDUCTION**

Using fluoroscopy, verify screw trajectories and confirm plate placement.

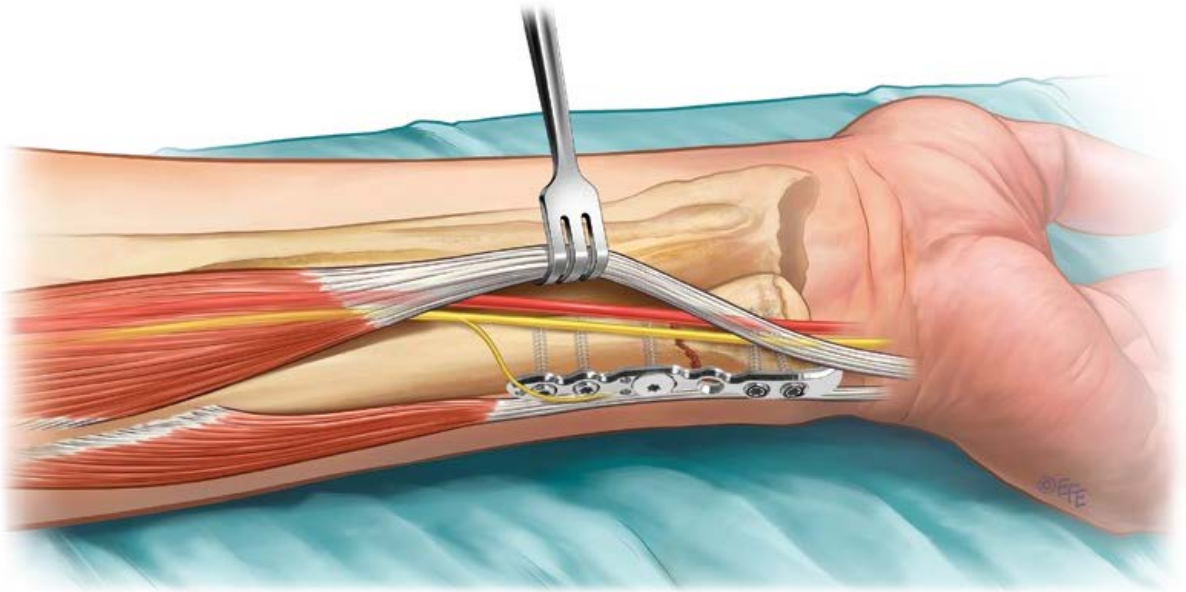

**Plate placement**



**Final placement confirmed**

## FINAL CONSTRUCT



Ulna Plates

# INSTRUMENT OVERVIEW

## RETRACTORS



Stabilizing Radiolucent Weitlaner 2x3, 5", Blunt Tip 6171.0004



Stabilizing Radiolucent Weitlaner 3x4, 8", Blunt Tip 6171.0005



Malleable Wire Replacement 6171.7008



Distal Radius Hohmann Retractor, 15mm 6171.7017



Senn Retractor 6171.7006



## ELEVATORS



Periosteal Elevator, Angled, Round Tip, 13mm 6171.7002

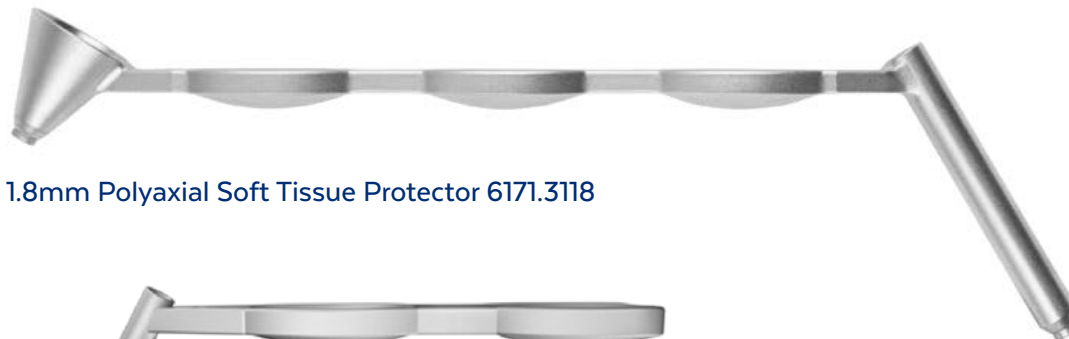


Distal Radius Freer Elevator 6171.7005

## DRILL GUIDES



1.8mm Speed Lock Drill Guide 6171.3218



1.8mm Polyaxial Soft Tissue Protector 6171.3118



Targeting Guide Drill Sleeve 6171.3512



1.8mm Threaded Drill Guide 6185.3218

## FORCEPS

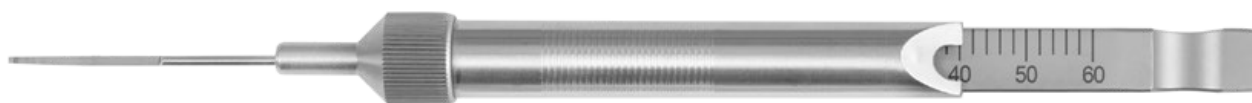


Point-to-Point Reduction Forceps,  
Ball Spikes, Ratcheting 6171.2003



Lobster Claw Reduction Forceps,  
Ratcheting 6179.2001

## DEPTH GAUGE



Depth Gauge, 60mm 6179.7020

## QUICK-CONNECT HANDLE



Small Handle, AO Quick-Connect 6171.7000

## DRIVER



T8 Driver, SR, 100mm, AO Quick-Connect 6179.6008

## K-WIRES



1.6mm K-wire, Trocar Tip, 150mm 6179.1116



1.6mm K-wire, Threaded Trocar Tip, 150mm 6171.1316

## TWIST DRILLS



1.8mm Drill Bit, 100mm, AO Quick-Connect 6171.5018



1.8mm Drill Bit, 130mm, AO Quick-Connect 6171.5019



1.5mm Drill Bit, 100mm, AO Quick-Connect 6171.5020



Calibrated 1.8mm Drill Bit, AO Quick-Connect, 100mm 6171.5518



Calibrated 1.8mm Drill Bit, AO Quick-Connect, 130mm 6171.5519

## HOOKS




Dental Pick, Curved Tip, Short Handle 6179.7012



Plate Positioning Tool 6171.7007
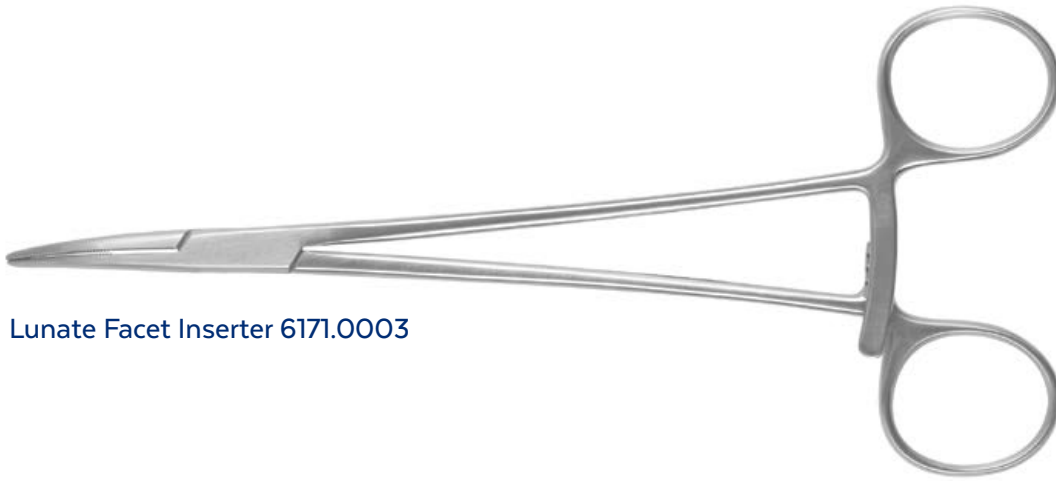
## LUNATE FACET HOOK INSTRUMENTS



Lunate Facet Inserter 6171.0003



Tamp 6171.7003




Mallet 603.977

## SCREW INSTRUMENTS



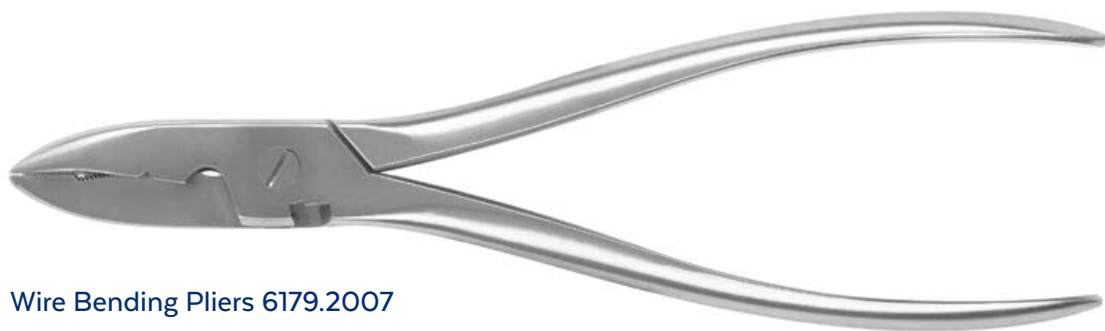Torque-Limiting Attachment, 1.2Nm, AO Quick-Connect 6171.5012

## PLIERS



Wire Bending Pliers 6179.2007

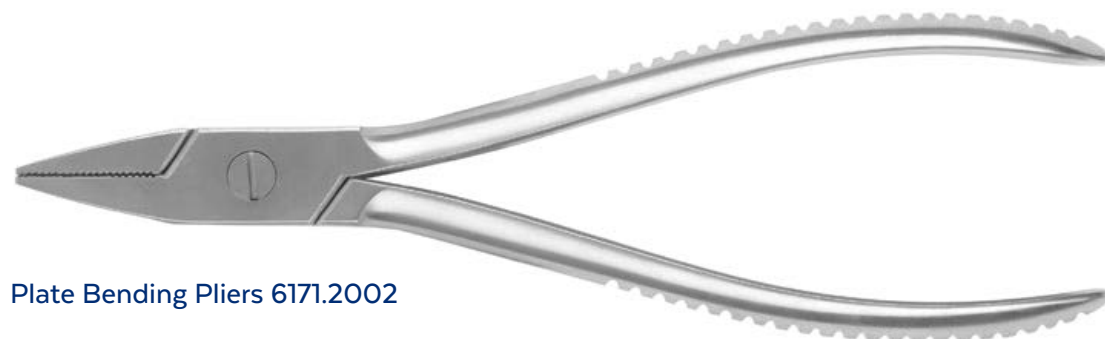


Plate Bending Pliers 6171.2002

## TARGETING GUIDES - ANTHEM® DOUBLE ROW VOLAR PLATES



Volar Plate Targeting Guide, Narrow, Left, 6171.3506

Volar Plate Targeting Guide, Standard, Left, 6171.3507

Volar Plate Targeting Guide, Wide, Left, 6171.3508



Volar Plate Targeting Guide, Narrow, Right, 6171.3509

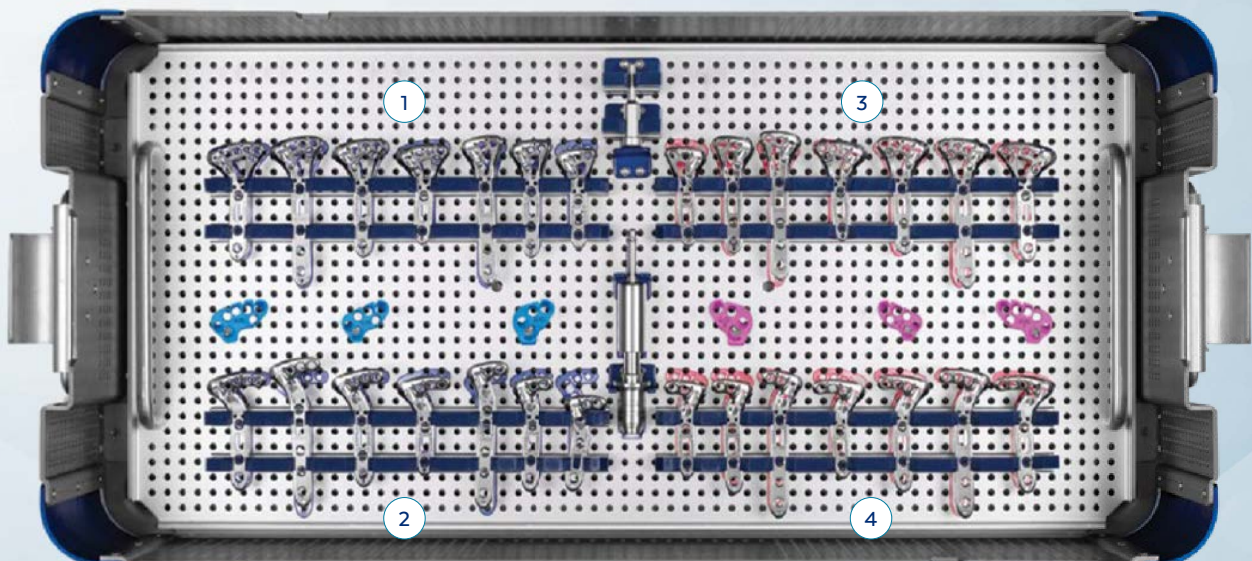
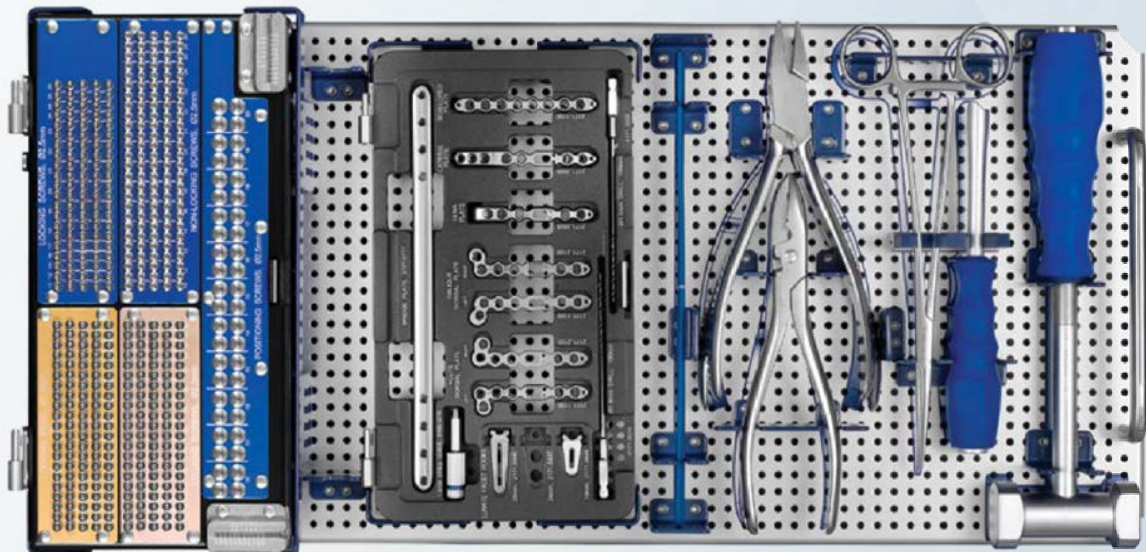
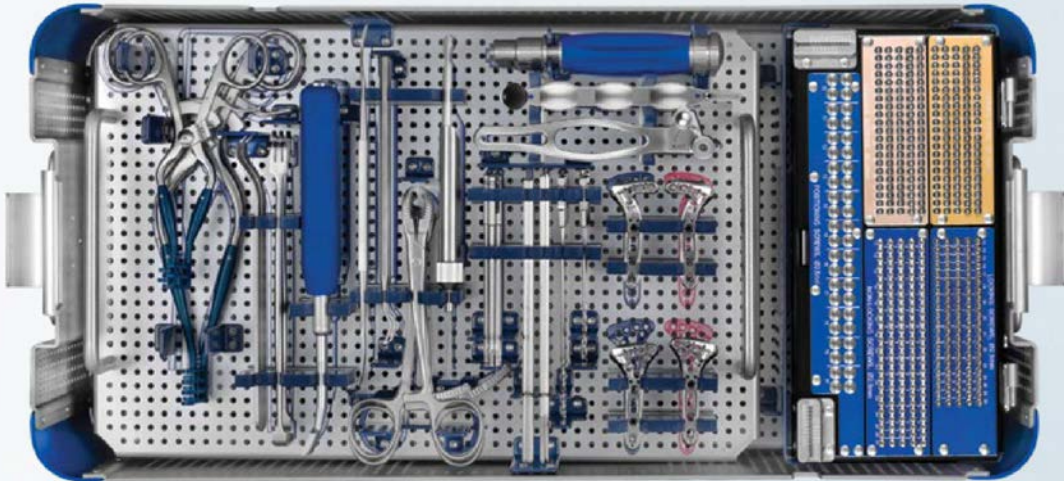
Volar Plate Targeting Guide, Standard, Right, 6171.3510

Volar Plate Targeting Guide, Wide, Right, 6171.3511

# ANTHEM<sup>®</sup> SS Distal Radius IMPLANTS 9171.9001

Part No.	Description	Qty
1	2171.1324 ANTHEM <sup>®</sup> Volar Plate, Double Row, Wide, Left, 4 Hole, SS	1
	2171.1313 ANTHEM <sup>®</sup> Volar Plate, Double Row, Standard, Left, 3 Hole, SS	1
	2171.1314 ANTHEM <sup>®</sup> Volar Plate, Double Row, Standard, Left, 4 Hole, SS	3
	2171.1315 ANTHEM <sup>®</sup> Volar Plate, Double Row, Standard, Left, 5 Hole, SS	1
	2171.1303 ANTHEM <sup>®</sup> Volar Plate, Double Row, Narrow, Left, 3 Hole, SS	1
	2171.1304 ANTHEM <sup>®</sup> Volar Plate, Double Row, Narrow, Left, 4 Hole, SS	1
	2171.1305 ANTHEM <sup>®</sup> Volar Plate, Double Row, Narrow, Left, 5 Hole, SS	1
2	2171.1024 ANTHEM <sup>®</sup> 7 Plate, Wide, Left, 4 Hole, SS	1
	2171.1013 ANTHEM <sup>®</sup> 7 Plate, Standard, Left, 3 Hole, SS	1
	2171.1014 ANTHEM <sup>®</sup> 7 Plate, Standard, Left, 4 Hole, SS	3
	2171.1015 ANTHEM <sup>®</sup> 7 Plate, Standard, Left, 5 Hole, SS	1
	2171.1003 ANTHEM <sup>®</sup> 7 Plate, Narrow, Left, 3 Hole, SS	1
	2171.1004 ANTHEM <sup>®</sup> 7 Plate, Narrow, Left, 4 Hole, SS	1
	2171.1005 ANTHEM <sup>®</sup> 7 Plate, Narrow, Left, 5 Hole, SS	1
3	2171.2303 ANTHEM <sup>®</sup> Volar Plate, Double Row, Narrow, Right, 3 Hole, SS	1
	2171.2304 ANTHEM <sup>®</sup> Volar Plate, Double Row, Narrow, Right, 4 Hole, SS	1
	2171.2305 ANTHEM <sup>®</sup> Volar Plate, Double Row, Narrow, Right, 5 Hole, SS	1
	2171.2313 ANTHEM <sup>®</sup> Volar Plate, Double Row, Standard, Right, 3 Hole, SS	1
	2171.2314 ANTHEM <sup>®</sup> Volar Plate, Double Row, Standard, Right, 4 Hole, SS	3
	2171.2315 ANTHEM <sup>®</sup> Volar Plate, Double Row, Standard, Right, 5 Hole, SS	1
	2171.2324 ANTHEM <sup>®</sup> Volar Plate, Double Row, Wide, Right, 4 Hole, SS	1
4	2171.2003 ANTHEM <sup>®</sup> 7 Plate, Narrow, Right, 3 Hole, SS	1
	2171.2004 ANTHEM <sup>®</sup> 7 Plate, Narrow, Right, 4 Hole, SS	1
	2171.2005 ANTHEM <sup>®</sup> 7 Plate, Narrow, Right, 5 Hole, SS	1
	2171.2013 ANTHEM <sup>®</sup> 7 Plate, Standard, Right, 3 Hole, SS	1
	2171.2014 ANTHEM <sup>®</sup> 7 Plate, Standard, Right, 4 Hole, SS	3
	2171.2015 ANTHEM <sup>®</sup> 7 Plate, Standard, Right, 5 Hole, SS	1
	2171.2024 ANTHEM <sup>®</sup> 7 Plate, Wide, Right, 4 Hole, SS	1
9171.0001	ANTHEM <sup>®</sup> SS Distal Radius System Graphic Case	










# ANTHEM® SS Distal Radius IMPLANTS 9171.9001 (CONT'D)

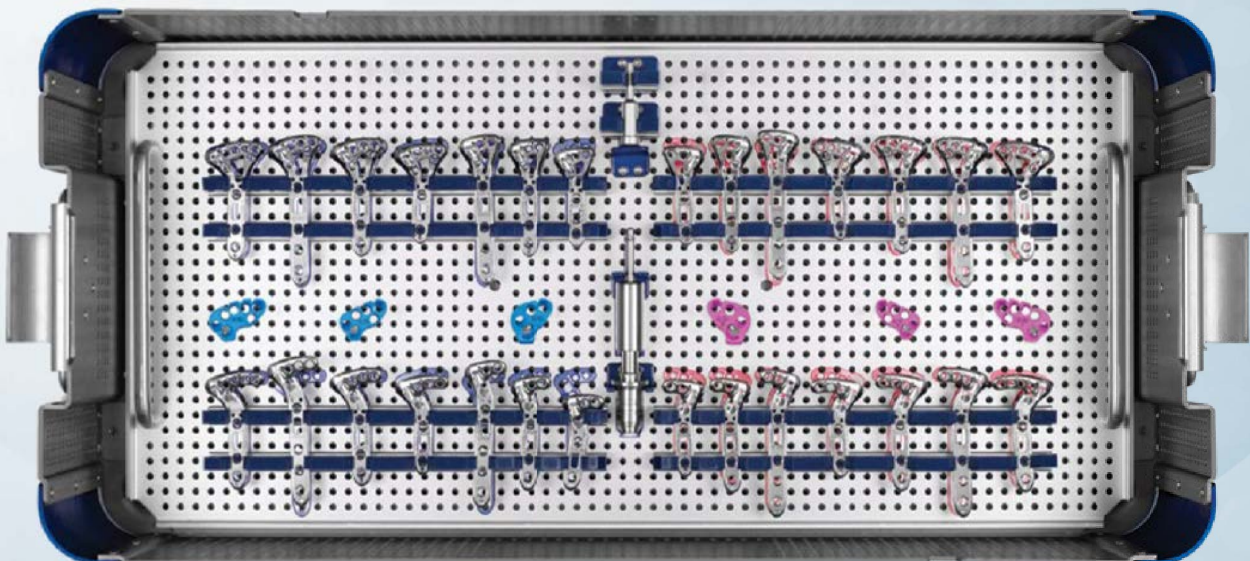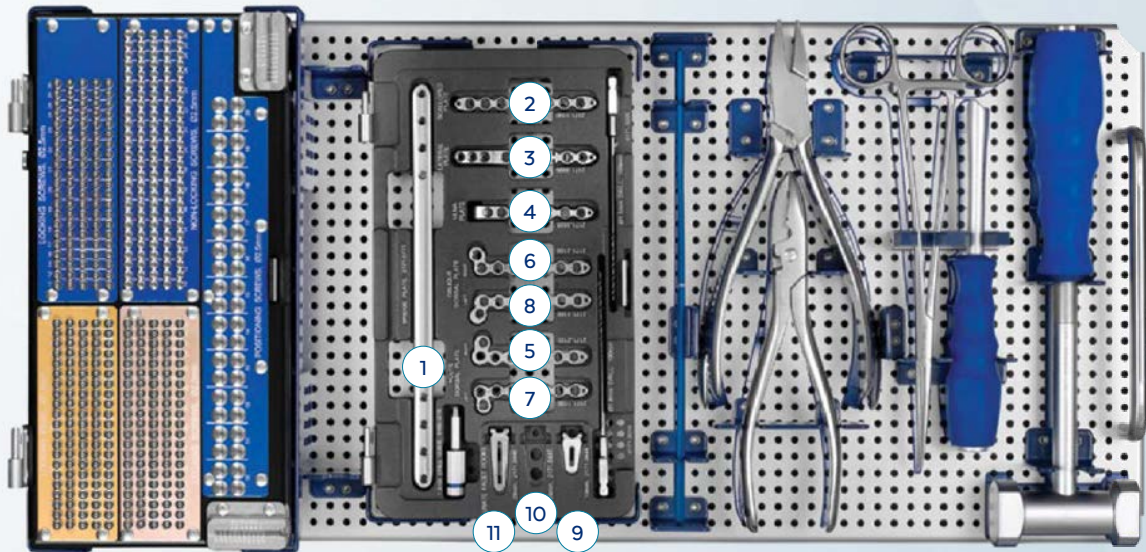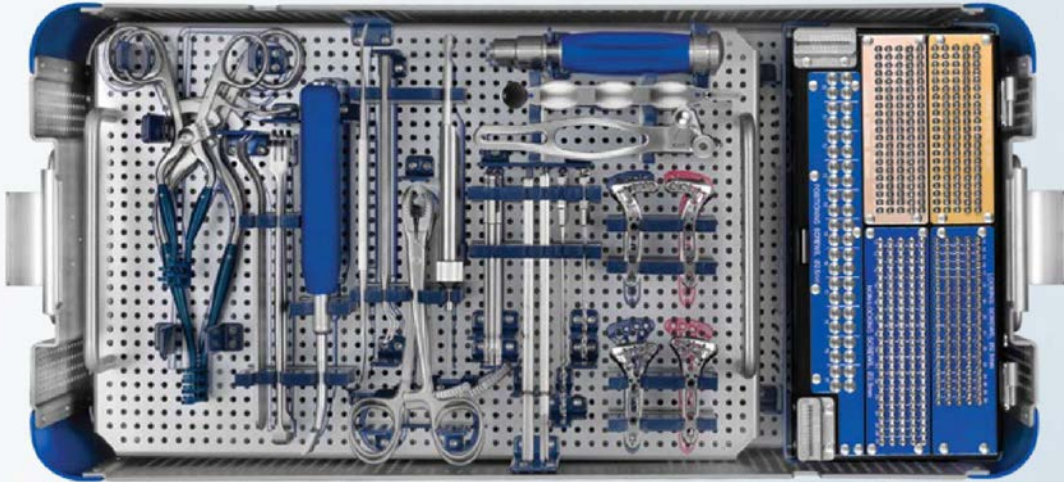	Part No.	Description	Qty
1	2171.0777	ANTHEM® Bridge Plate, SS	1
2	2171.1100	ANTHEM® Scalloped Plate, SS	1
3	2171.0999	ANTHEM® Lateral Plate, SS	1
4	2171.0555	ANTHEM® Ulna Plate, SS	1
5	2171.2155	ANTHEM® Dorsal Plate, Right, Acute, SS	1
6	2171.2188	ANTHEM® Dorsal Plate, Right, Oblique, SS	1
7	2171.1155	ANTHEM® Dorsal Plate, Left, Acute, SS	1
8	2171.1188	ANTHEM® Dorsal Plate, Left, Oblique, SS	1
9	2171.0444	ANTHEM® Lunate Facet Hook, 19mm, SS	1
10	2171.0045	ANTHEM® Lunate Facet Hook, 24mm, SS	1
11	2171.0446	ANTHEM® Lunate Facet Hook, 29mm, SS	1
	9171.1001	ANTHEM® SS Distal Radius System Fragment Specific Module	







# ANTHEM® Ti Distal Radius IMPLANTS 9171.9002

Part No.	Description	Qty
1171.1324	ANTHEM® Volar Plate, Double Row, Wide, Left, 4 Hole, Ti	1
1171.1313	ANTHEM® Volar Plate, Double Row, Standard, Left, 3 Hole, Ti	1
1171.1314	ANTHEM® Volar Plate, Double Row, Standard, Left, 4 Hole, Ti	3
1171.1315	ANTHEM® Volar Plate, Double Row, Standard, Left, 5 Hole, Ti	1
1171.1303	ANTHEM® Volar Plate, Double Row, Narrow, Left, 3 Hole, Ti	1
1171.1304	ANTHEM® Volar Plate, Double Row, Narrow, Left, 4 Hole, Ti	1
1171.1305	ANTHEM® Volar Plate, Double Row, Narrow, Left, 5 Hole, Ti	1
1171.1024	ANTHEM® 7 Plate, Wide, Left, 4 Hole, Ti	1
1171.1013	ANTHEM® 7 Plate, Standard, Left, 3 Hole, Ti	1
1171.1014	ANTHEM® 7 Plate, Standard, Left, 4 Hole, Ti	3
1171.1015	ANTHEM® 7 Plate, Standard, Left, 5 Hole, Ti	1
1171.1003	ANTHEM® 7 Plate, Narrow, Left, 3 Hole, Ti	1
1171.1004	ANTHEM® 7 Plate, Narrow, Left, 4 Hole, Ti	1
1171.1005	ANTHEM® 7 Plate, Narrow, Left, 5 Hole, Ti	1
1171.2303	ANTHEM® Volar Plate, Double Row, Narrow, Right, 3 Hole, Ti	1
1171.2304	ANTHEM® Volar Plate, Double Row, Narrow, Right, 4 Hole, Ti	1
1171.2305	ANTHEM® Volar Plate, Double Row, Narrow, Right, 5 Hole, Ti	1
1171.2313	ANTHEM® Volar Plate, Double Row, Standard, Right, 3 Hole, Ti	1
1171.2314	ANTHEM® Volar Plate, Double Row, Standard, Right, 4 Hole, Ti	3
1171.2315	ANTHEM® Volar Plate, Double Row, Standard, Right, 5 Hole, Ti	1
1171.2324	ANTHEM® Volar Plate, Double Row, Wide, Right, 4 Hole, Ti	1
1171.2003	ANTHEM® 7 Plate, Narrow, Right, 3 Hole, Ti	1
1171.2004	ANTHEM® 7 Plate, Narrow, Right, 4 Hole, Ti	1
1171.2005	ANTHEM® 7 Plate, Narrow, Left, 5 Hole, Ti	1
1171.2013	ANTHEM® 7 Plate, Standard, Right, 3 Hole, Ti	1
1171.2014	ANTHEM® 7 Plate, Standard, Right, 4 Hole, Ti	3
1171.2015	ANTHEM® 7 Plate, Standard, Right, 5 Hole, Ti	1
1171.2024	ANTHEM® 7 Plate, Wide, Right, 4 Hole, Ti	1
9171.0002	ANTHEM® Ti Distal Radius System Graphic Case	

# ANTHEM® Ti Distal Radius IMPLANTS 9171.9002 (CONT'D)

Part No.	Description	Qty
1171.0777	ANTHEM® Bridge Plate, Ti	1
1171.1100	ANTHEM® Scalloped Plate, Ti	1
1171.0999	ANTHEM® Lateral Plate, Ti	1
1171.0555	ANTHEM® Ulna Plate, Ti	1
1171.2155	ANTHEM® Dorsal Plate, Right, Acute, Ti	1
1171.2188	ANTHEM® Dorsal Plate, Right, Oblique, Ti	1
1171.1155	ANTHEM® Dorsal Plate, Left, Acute, Ti	1
1171.1188	ANTHEM® Dorsal Plate, Left, Oblique, Ti	1
1171.0444	ANTHEM® Lunate Facet Hook, 19mm, Ti	1
1171.0445	ANTHEM® Lunate Facet Hook, 24mm, Ti	1
1171.0446	ANTHEM® Lunate Facet Hook, 29mm, Ti	1
9171.1002	ANTHEM® Ti Distal Radius System Fragment Specific Module	



# ANTHEM<sup>®</sup> Distal Radius IMPLANTS 9171.9003

## 1 1.8mm Locking Pegs

Part No.	Diameter/Length	Qty
7171.4810	1.8x10mm, CoCr	5
7171.4812	1.8x12mm, CoCr	5
7171.4814	1.8x14mm, CoCr	5
7171.4816	1.8x16mm, CoCr	5
7171.4817	1.8x17mm, CoCr	5
7171.4818	1.8x18mm, CoCr	8
7171.4819	1.8x19mm, CoCr	8
7171.4820	1.8x20mm, CoCr	8
7171.4821	1.8x21mm, CoCr	8
7171.4822	1.8x22mm, CoCr	8
7171.4823	1.8x23mm, CoCr	8
7171.4824	1.8x24mm, CoCr	4
7171.4826	1.8x26mm, CoCr	4
7171.4828	1.8x28mm, CoCr	4
7171.4830	1.8x30mm, CoCr	2

## 2 2.0mm Locking Screws

Part No.	Diameter/Length	Qty
7171.5010	2.0x10mm, CoCr	2
7171.5012	2.0x12mm, CoCr	2
7171.5014	2.0x14mm, CoCr	2
7171.5016	2.0x16mm, CoCr	2
7171.5017	2.0x17mm, CoCr	2
7171.5018	2.0x18mm, CoCr	8
7171.5019	2.0x19mm, CoCr	8
7171.5020	2.0x20mm, CoCr	8
7171.5021	2.0x21mm, CoCr	8
7171.5022	2.0x22mm, CoCr	8
7171.5023	2.0x23mm, CoCr	8
7171.5024	2.0x24mm, CoCr	4
7171.5026	2.0x26mm, CoCr	4
7171.5028	2.0x28mm, CoCr	4
7171.5030	2.0x30mm, CoCr	2

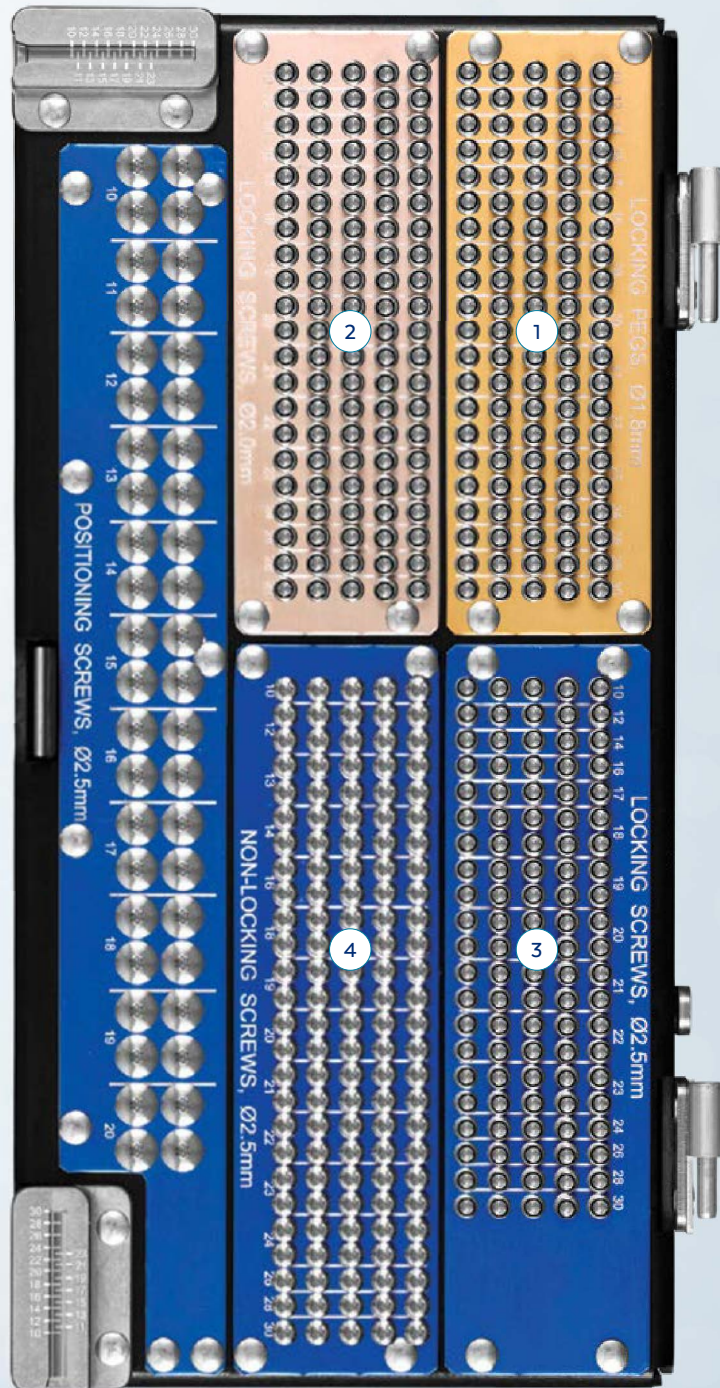
## 3 2.5mm Locking Screws

Part No.	Diameter/Length	Qty
7171.5510	2.5x10mm, CoCr	5
7171.5512	2.5x12mm, CoCr	5
7171.5514	2.5x14mm, CoCr	5
7171.5516	2.5x16mm, CoCr	5
7171.5517	2.5x17mm, CoCr	5
7171.5518	2.5x18mm, CoCr	8
7171.5519	2.5x19mm, CoCr	8
7171.5520	2.5x20mm, CoCr	8
7171.5521	2.5x21mm, CoCr	8
7171.5522	2.5x22mm, CoCr	8
7171.5523	2.5x23mm, CoCr	8
7171.5524	2.5x24mm, CoCr	4
7171.5526	2.5x26mm, CoCr	4
7171.5528	2.5x28mm, CoCr	4
7171.5530	2.5x30mm, CoCr	2

## 4 2.5mm Non-Locking Screws

Part No.	Diameter/Length	Qty
7171.6510	2.5x10mm, CoCr	2
7171.6512	2.5x12mm, CoCr	10
7171.6514	2.5x14mm, CoCr	10
7171.6516	2.5x16mm, CoCr	10
7171.6518	2.5x18mm, CoCr	10
7171.6519	2.5x19mm, CoCr	10
7171.6520	2.5x20mm, CoCr	10
7171.6521	2.5x21mm, CoCr	10
7171.6522	2.5x22mm, CoCr	10
7171.6523	2.5x23mm, CoCr	10
7171.6524	2.5x24mm, CoCr	10
7171.6526	2.5x26mm, CoCr	2
7171.6528	2.5x28mm, CoCr	2
7171.6530	2.5x30mm, CoCr	2






# ANTHEM<sup>®</sup> Distal Radius IMPLANTS 9171.9003 (CONT'D)

## 1 2.5mm Positioning Screws

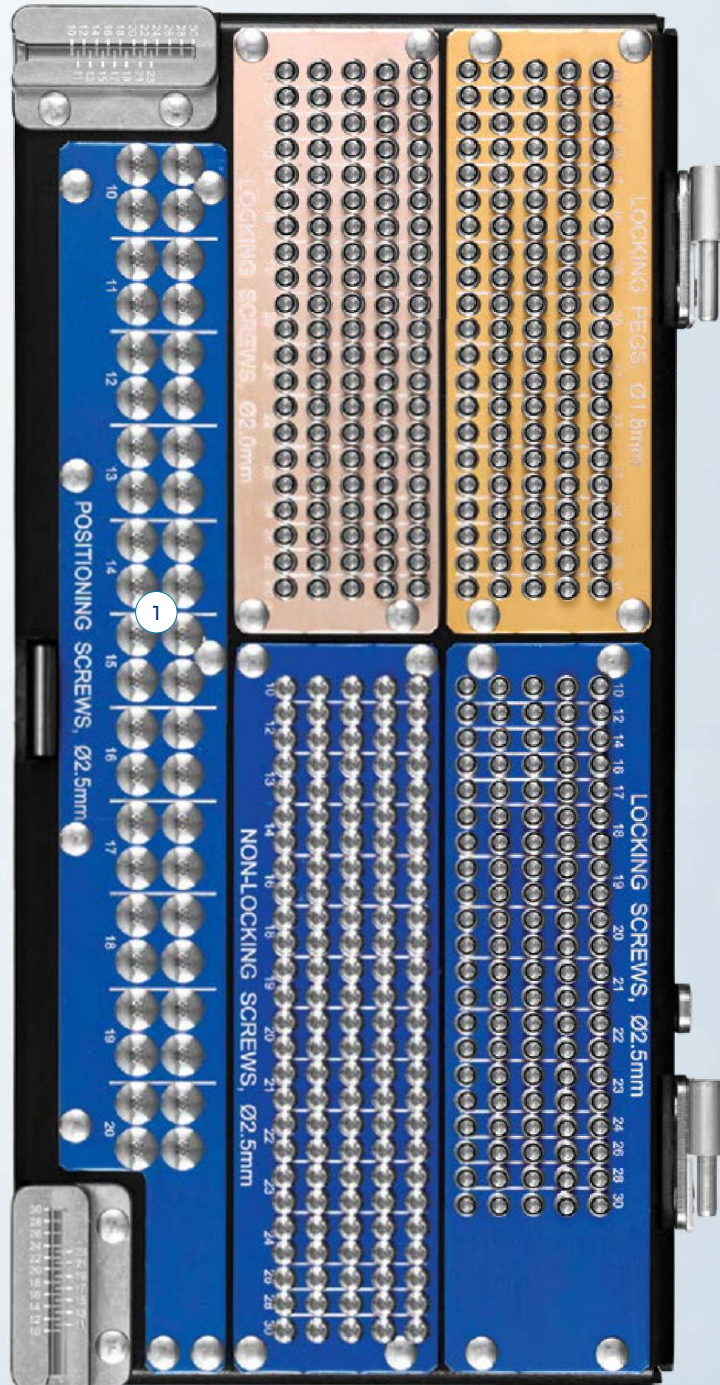Part No.	Diameter/Length	Qty
7171.7510	2.5x10mm, CoCr	4
7171.7511	2.5x11mm, CoCr	4
7171.7512	2.5x12mm, CoCr	4
7171.7513	2.5x13mm, CoCr	4
7171.7514	2.5x14mm, CoCr	4
7171.7515	2.5x15mm, CoCr	4
7171.7516	2.5x16mm, CoCr	4
7171.7517	2.5x17mm, CoCr	4
7171.7518	2.5x18mm, CoCr	4
7171.7519	2.5x19mm, CoCr	4
7171.7520	2.5x20mm, CoCr	4
9171.9003	ANTHEM <sup>®</sup> Distal Radius Screw Module	





# ANTHEM® Distal Radius Fracture System

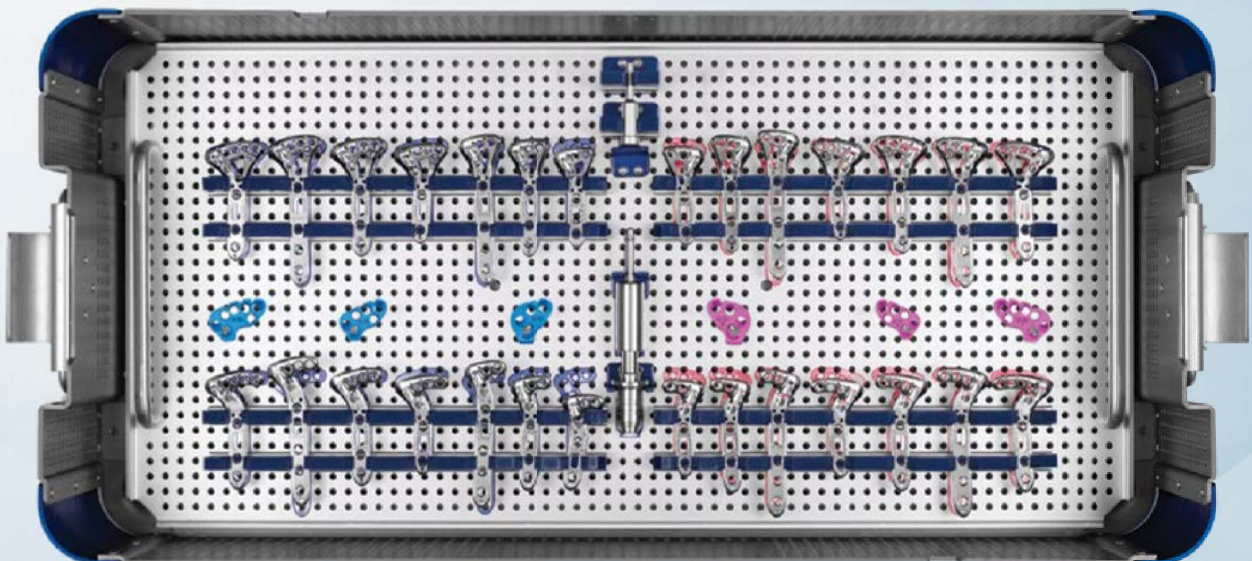
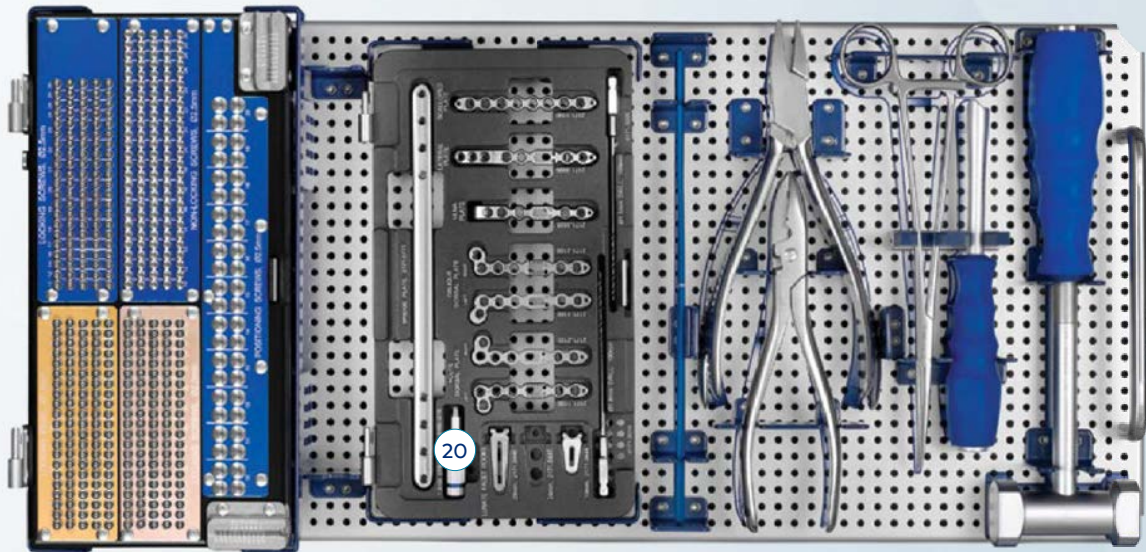
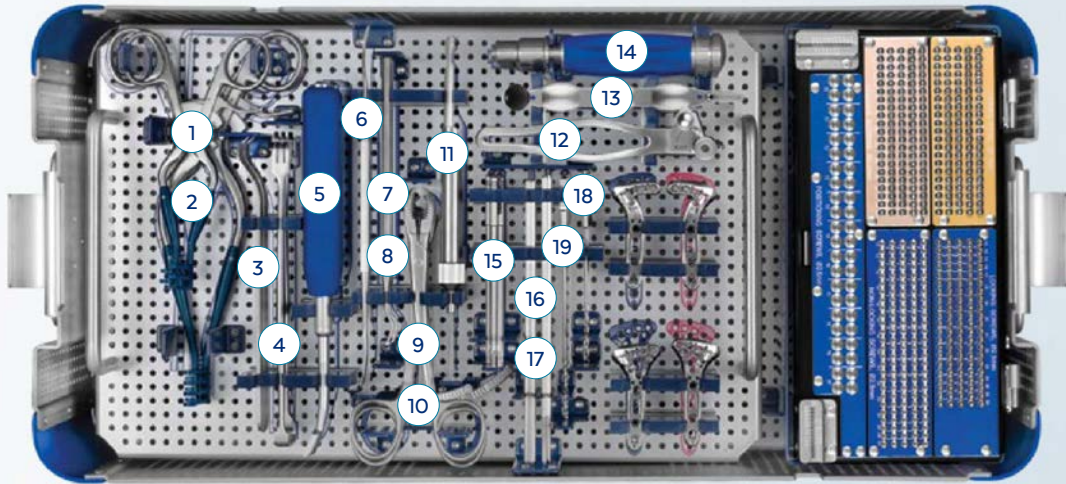
## INSTRUMENTS 9171.9001 AND 9171.9002

	Part No.	Description	Qty
1	6171.0004	Stabilizing Radiolucent Weitlaner 2x3 5", Blunt Tip	1
2	6171.0005	Stabilizing Radiolucent Weitlaner 3x4 8", Blunt Tip	1
3	6171.7017	Distal Radius Hohmann Retractor, 15mm	2
4	6171.7006	Senn Retractor	2
5	6171.7002	Periosteal Elevator, Angled, Round Tip, 13mm	1
6	6171.7005	Distal Radius Freer Elevator	1
7	6179.7012	Dental Pick, Curved Tip, Short Handle	1
8	6171.7007	Plate Positioning Tool	1
9	6179.2001	Lobster Claw Reduction Forceps, Ratcheting	2
10	6171.2003	Point-to-Point Reduction Forceps, Ball Spikes, Ratcheting	2
11	6179.7020	Depth Gauge, 60mm	1
12	6171.3218	1.8mm Speed Lock Drill Guide	1
13	6171.3118	1.8mm Polyaxial Soft Tissue Protector	1
14	6171.7000	Small Handle, AO Quick-Connect	2
15	6179.6008	T8 Driver, SR, 100mm, AO Quick-Connect	5
16	6179.1116	1.6mm K-Wire, Trocar Tip, 150mm	10
17	6171.1316	1.6mm K-Wire, Threaded Trocar Tip, 150mm	5
18	6171.5018	1.8mm Drill Bit, 100mm, AO Quick-Connect	5
19	6171.5019	1.8mm Drill Bit, 130mm, AO Quick-Connect	5
20	6185.3218	1.8mm Threaded Drill Guide	2

### Additionally Available

6171.3110	1.8mm Small Polyaxial Drill Guide
6171.7009	Small Torque-Limiting Handle, 1.2Nm, AO Quick-Connect










# ANTHEM<sup>®</sup> Distal Radius Fracture System

## INSTRUMENTS 9171.9001 AND 9171.9002 (CONT'D)

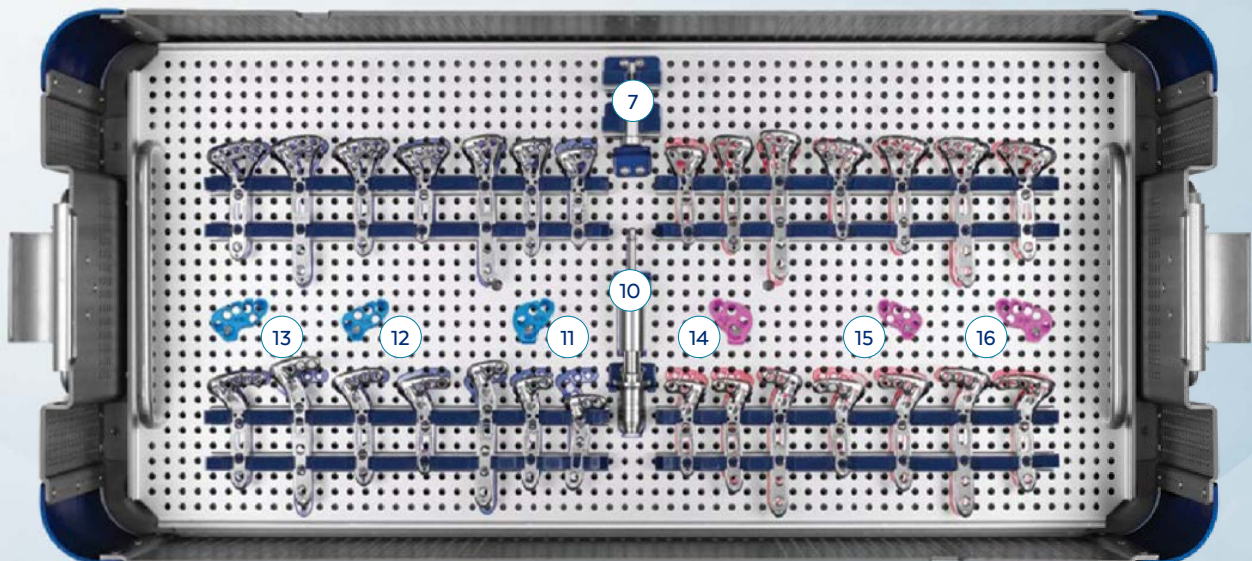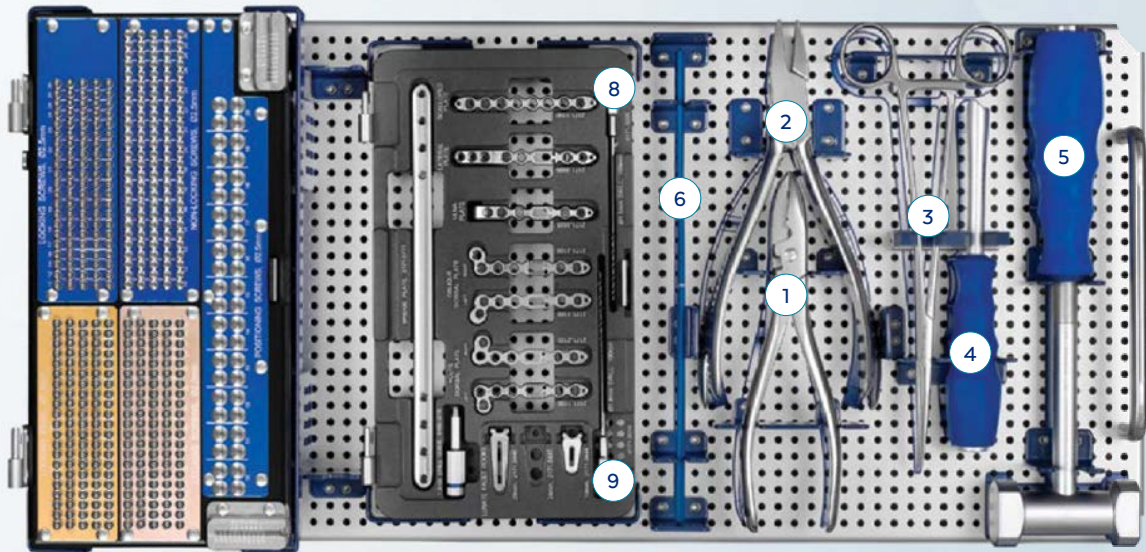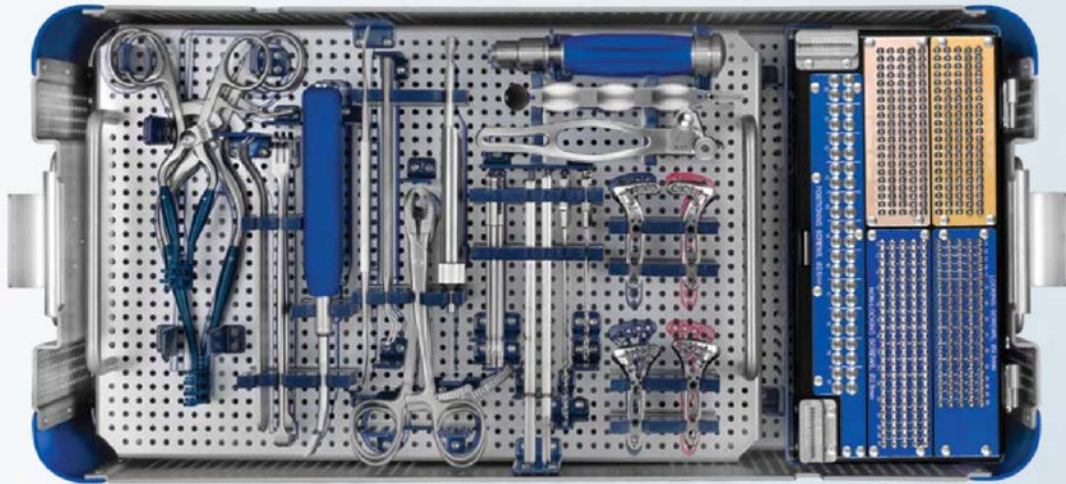	Part No.	Description	Qty
1	6179.2007	Wire Bending Pliers	1
2	6171.2002	Plate Bending Pliers	2
3	6171.0003	Lunate Facet Inserter	1
4	6171.7003	Tamp	1
5	603.977	Mallet	1
6	6171.7008	Malleable Wire Replacement	5
7	6171.3512	Targeting Guide Drill Sleeve	1
8	6171.5020	1.5mm Drill Bit, 100mm, AO Quick-Connect	2
9	6171.5018	1.8mm Drill Bit, 100mm, AO Quick-Connect	2
10	6171.5012	Torque-Limiting Attachment, 1.2Nm, AO Quick-Connect	1
11	6171.3506	Volar Plate Targeting Guide, Narrow, Left	1
12	6171.3507	Volar Plate Targeting Guide, Standard, Left	1
13	6171.3508	Volar Plate Targeting Guide, Wide, Left	1
14	6171.3509	Volar Plate Targeting Guide, Narrow, Right	1
15	6171.3510	Volar Plate Targeting Guide, Standard, Right	1
16	6171.3511	Volar Plate Targeting Guide, Wide, Right	1
	6171.3513	Targeting Guide Depth Gauge Cap	1







# IMPORTANT INFORMATION ON ANTHEM® FRACTURE SYSTEM

## DESCRIPTION

The ANTHEM® Fracture System is a family of plates and screws designed to be used for internal bone fixation. The implants are available in various sizes and shapes to accommodate patient anatomy, and may be contoured or straight, with locking and non-locking screws. ANTHEM® implants are manufactured from titanium, titanium alloy, cobalt chromium molybdenum alloy, or stainless steel, as specified in ASTM F67, F136, F1295, F1472, F1537, F2229, F138 and F139. All implants are for single use only.

## INDICATIONS

The ANTHEM® Fracture System is indicated for fixation of fractures, osteotomies, arthrodesis and reconstruction of bones for the appropriate size of the device to be used in adult patients, including the clavicle, scapula, humerus, radius, ulna, small bones (metacarpals, metatarsals, phalanges), wrist, pelvis, femur, tibia, fibula, ankle, and foot. The clavicle hook plate may be used for dislocations of the acromioclavicular joint. Mini fragment plates are also indicated for fixation of fractures of the acetabulum, patella, and bone fragments, replantation, malunions and nonunion, and for non-load bearing stabilization and reduction of long bone fragments. Metaphyseal plates are indicated for non-load bearing stabilization and reduction of long bone fragments, and for fixation of bones including the radius and ulna.

Small fragment, mini fragment, proximal tibia, clavicle, metaphyseal, and distal fibula plates may be used in all pediatric subgroups (except neonates) and small stature adults. Distal radius, distal tibia, metaphyseal, and mini fragment plates may be used in adolescents (12-21 years of age). Plating may be used in patients with osteopenic bone.

## CONTRAINDICATIONS

Use of these implants is contraindicated in patients with the following conditions:

- Any active or suspended latent infection or marked local inflammation in or about the affected area.
- Compromised vascularity that would inhibit adequate blood supply to the fracture or the operative site.
- Bone stock compromised by disease, infection or prior implantation that cannot provide adequate support and/or fixation of the devices.
- Use of plating on or around growth plates in pediatric patients.
- Material sensitivity, documented or suspected.
- Obesity. An overweight or obese patient can produce loads on the implant that can lead to failure of the device itself.
- Patients having inadequate tissue coverage over the operative site.
- Implant utilization that would interfere with anatomical structures or physiological performance.
- Any mental or neuromuscular disorder which would create an unacceptable risk of fixation failure or complications in postoperative care.
- Other medical or surgical conditions which would preclude the potential benefit of surgery.

## WARNINGS

The correct implant selection is extremely important. Failure to use the appropriate implant for the fracture condition may accelerate clinical failure. Failure to use the proper component to maintain adequate blood supply and provide rigid fixation may result in loosening, bending, cracking or fracture of the implant and/or bone. The correct implant size for a given patient can be determined by evaluating the patient's height, weight, functional demands and anatomy. Every implant must be used in the correct anatomic location, consistent with accepted standards of internal fixation.

## PRECAUTIONS

The implantation of fixation devices should be performed only by experienced surgeons with specific training in the use of this system because this is a technically demanding procedure presenting a risk of serious injury to the patient. Preoperative planning and patient anatomy should be considered when selecting implant size.

Surgical implants must never be reused. Even though the device appears undamaged, it may have small defects and internal stress patterns which could lead to breakage.

## MR SAFETY INFORMATION

These devices have not been evaluated for safety and compatibility in the MR environment. It has not been tested for heating, migration, or image artifact in the MR environment. The safety of these devices in the MR environment is unknown. Scanning a patient who has this device may result in patient injury.

## CAUTIONS

### Pre-operative

- These implants are for single use only.

- Implants that came in contact with body fluids should never be reused.
- Ensure that all components needed for surgery are available in the surgical suite.
- Inspection is recommended prior to surgery to determine if implants have been damaged during storage.
- While rare, intra-operative fracture or breakage of instruments can occur. Instruments which have experienced excessive use or excessive force are susceptible to fracture. Instruments should be examined for wear or damage prior to surgery.

### Intra-operative

- Avoid surface damage of implants.
- Discard all damaged or mishandled implants.
- Contouring or bending of an implant should be avoided where possible, because it may reduce its fatigue strength and can cause failure under load.
- Implants are available in different versions, varying for example in length, diameter, material and number of drilled holes. Select the required version carefully.
- During the course of the operation, repeatedly check to ensure that the connection between the implant and the instrument, or between the instruments, is secure.
- Implants which consist of several components must only be used in the prescribed combination (refer to the ANTHEM® Surgical Technique Guide).
- After the procedure check the proper positioning of all implants using the image intensifier.
- Do not use components from this system in conjunction with components from any other manufacturer's system unless otherwise specified (refer to the ANTHEM® Surgical Technique Guide).

### Post-operative

- Post-operative patient activity: These implants are neither intended to carry the full load of the patient acutely, nor intended to carry a significant portion of the load for extended periods of time. For this reason post-operative instructions and warnings to patients are extremely important. External immobilization (e.g. bracing or casting) may be employed until X-rays or other procedures confirm adequate bone consolidation.
- The implant is a short-term implant. In the event of a delay in bone consolidation, or if such consolidation does not take place, or if explantation is not carried out, complications may occur, for example fracture or loosening of the implant or instability of the implant system. Regular post-operative examinations (e.g., X-ray checks) are advisable.
- The risk of post-operative complication (e.g. failure of an implant) is higher if patients are obese and/or cannot follow the recommendations of the physician because of any mental or neuromuscular disorder. For this reason those patients must have additional post-operative follow-up.
- Implant removal should be followed by adequate postoperative management to avoid fracture or refracture of the bone.

### Informing the Patient

The implant affects the patient's ability to carry loads and her/his mobility and general living circumstances. The surgeon must counsel each patient individually on correct behavior and activity after the implantation.

The surgeon must warn each patient that the device cannot and does not replicate a normally healthy bone, that the device can break or become damaged as a result of strenuous activity, trauma, mal-union or non-union and that the device has a finite expected service life and may need to be removed at some time in the future.

## ADVERSE EFFECTS

In many instances, adverse results may be clinically related rather than device related. The following are the most frequent adverse effects involving the use of internal fracture fixation devices:

- Delayed union or non-union of the fracture site.
- These devices can break when subjected to the increased loading associated with delayed unions and/or non-unions. Internal fixation devices are load sharing devices which are intended to hold fracture bone surface in a position to facilitate healing. If healing is delayed or does not occur, the appliance may eventually break due to metal fatigue. Loads on the device produced by load bearing and the patient's activity level will dictate the longevity of the device.
- Conditions attributable to non-union, osteoporosis, osteomalacia, diabetes, inhibited revascularization and poor bone formation can cause loosening, bending, cracking, fracture of the device or premature loss of rigid fixation with the bone.

## IMPORTANT INFORMATION ON ANTHEM® FRACTURE SYSTEM

- Improper alignment can cause a mal-union of the bone and/or bending, cracking or even breakage of the device.
- Increased fibrous tissue response around the fracture site due to unstable comminuted fractures.
- Early or late infection, deep or superficial.
- Deep venous thrombosis.
- Avascular necrosis.
- Shortening of the effected bone/fracture site.
- Subclinical nerve damage may possibly occur as a result of the surgical trauma.
- Material sensitivity reactions in patients following surgical implantation have rarely been reported, however their significance awaits further clinical evaluation.

### PACKAGING

These implants may be supplied pre-packaged and sterile, using gamma irradiation. The integrity of the sterile packaging should be checked to ensure that sterility of the contents is not compromised. Packaging should be carefully checked for completeness and all components should be carefully checked to ensure that there is no damage prior to use. Damaged packages or products should not be used, and should be returned to Globus Medical. During surgery, after the correct size has been determined, remove the products from the packaging using aseptic technique.

The instruments are provided nonsterile and are steam sterilized prior to use, as described in the STERILIZATION section below. Following use or exposure to soil, instruments and instrument trays and cases must be cleaned, as described in the CLEANING section below.

### HANDLING

All instruments and implants should be treated with care. Improper use or handling may lead to damage and/or possible malfunction. Instruments should be checked to ensure that they are in working order prior to surgery.

Implants are single use devices and should not be cleaned. Re-cleaning of single use implants might lead to mechanical failure and/or material degradation. Discard any implants that may have been accidentally contaminated.

### CLEANING

Instruments should be cleaned separately from instrument trays and cases. Lids should be removed from cases for the cleaning process, if applicable. All instruments that can be disassembled must be disassembled for cleaning. All handles must be detached. Instruments may be reassembled following sterilization. The products should be cleaned using neutral cleaners before sterilization and introduction into a sterile surgical field or (if applicable) return of the product to Globus Medical.

Cleaning and disinfecting can be performed with aldehyde-free solvents at higher temperatures. Cleaning and decontamination must include the use of neutral cleaners followed by a deionized water rinse. Note: certain cleaning solutions such as those containing formalin, glutaraldehyde, bleach and/or other alkaline cleaners may damage some devices, particularly instruments; these solutions should not be used.

The following cleaning methods should be observed when cleaning instruments and instrument trays and cases after use or exposure to soil, and prior to sterilization:

1. Immediately following use, ensure that the instruments are wiped down to remove all visible soil and kept from drying by submerging or covering with a wet towel.
2. Disassemble all instruments that can be disassembled.
3. Rinse the instruments under running tap water to remove all visible soil. Flush the lumens a minimum of 3 times, until the lumens flush clean.
4. Prepare Enzo® (or a similar enzymatic detergent) per manufacturer's recommendations.
5. Immerse the instruments in the detergent and allow them to soak for a minimum of 2 minutes.
6. Use a soft bristled brush to thoroughly clean the instruments. Use a pipe cleaner for any lumens. Pay close attention to hard to reach areas.
7. Using a sterile syringe, draw up the enzymatic detergent solution. Flush any lumens and hard to reach areas until no soil is seen exiting the area.
8. Remove the instruments from the detergent and rinse them in running warm tap water.
9. Prepare Enzo® (or a similar enzymatic detergent) per manufacturer's recommendations in an ultrasonic cleaner.

10. Completely immerse the instruments in the ultrasonic cleaner and ensure detergent is in lumens by flushing the lumens. Sonicate for a minimum of 3 minutes.

11. Remove the instruments from the detergent and rinse them in running deionized water or reverse osmosis water for a minimum of 2 minutes.

12. Dry instruments using a clean soft cloth and filtered pressurized air.

13. Visually inspect each instrument for visible soil. If visible soil is present, then repeat cleaning process starting with Step 3.

### CONTACT INFORMATION

Globus Medical may be contacted at 1-866-GLOBUS1 (456-2871). A surgical technique manual may be obtained by contacting Globus Medical.

### STERILIZATION

These implants may be available sterile or nonsterile. Instruments are available nonsterile.

Sterile implants are sterilized by gamma radiation, validated to ensure a Sterility Assurance Level (SAL) of 10<sup>-6</sup>. Sterile products are packaged in a heat sealed, Tyvek pouch or in a container/pouch. The expiration date is provided in the package label. These products are considered sterile unless the packaging has been opened or damaged. Sterile implants meet pyrogen limit specifications.

Nonsterile implants and instruments have been validated to ensure an SAL of 10<sup>-6</sup>. The use of an FDA-cleared wrap is recommended, per the Association for the Advancement of Medical Instrumentation (AAMI) ST79, *Comprehensive Guide to Steam Sterilization and Sterility Assurance in Health Care Facilities*. It is the end user's responsibility to use only sterilizers and accessories (such as sterilization wraps, sterilization pouches, chemical indicators, biological indicators, and sterilization cassettes) that have been cleared by the FDA for the selected sterilization cycle specifications (time and temperature).

When using a rigid sterilization container, the following must be taken into consideration for proper sterilization of Globus devices and loaded graphic cases:


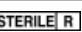


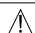

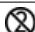



- Recommended sterilization parameters are listed in the table below.
- Only FDA-cleared rigid sterilization containers for use with pre-vacuum steam sterilization may be used.
- When selecting a rigid sterilization container, it must have a minimum filter area of 176 in<sup>2</sup> total, or a minimum of four (4) 7.5in diameter filters.
- No more than one (1) loaded graphic case or its contents can be placed directly into a rigid sterilization container.
- Stand-alone modules/racks or single devices must be placed, without stacking, in a container basket to ensure optimal ventilation.
- The rigid sterilization container manufacturer's instructions for use are to be followed; if questions arise, contact the manufacturer of the specific container for guidance.
- Refer to AAMI ST79 for additional information concerning the use of rigid sterilization containers.

For implants and instruments provided NONSTERILE, sterilization is recommended (wrapped or containerized) as follows:

Method	Cycle Type	Temperature	Exposure Time	Drying Time
Steam	Pre-vacuum	132°C (270°F)	4 Minutes	30 Minutes

*These parameters are validated to sterilize only this device. If other products are added to the sterilizer, the recommended parameters are not valid and new cycle parameters must be established by the user. The sterilizer must be properly installed, maintained, and calibrated. Ongoing testing must be performed to confirm inactivation of all forms of viable microorganisms.*

**CAUTION:** Federal (USA) Law Restricts this Device to Sale by or on the order of a Physician.

SYMBOL TRANSLATION			
	CATALOGUE NUMBER		STERILIZED BY IRRADIATION
	LOT NUMBER		AUTHORISED REPRESENTATIVE IN THE EUROPEAN COMMUNITY
	CAUTION		MANUFACTURER
	SINGLE USE ONLY		USE BY (YYYY-MM-DD)
	QUANTITY		PRESCRIPTION USE ONLY

DI201A REV E



Globus Medical  
Valley Forge Business Center  
2560 General Armistead Avenue  
Audubon, PA 19403  
[www.globusmedical.com](http://www.globusmedical.com)

Customer Service:

Phone 1-866-GLOBUS1 (or 1-866-456-2871)

Fax 1-866-GLOBUS3 (or 1-866-456-2873)

©2020 Globus Medical. All rights reserved. Patent [www.globusmedical.com/patents](http://www.globusmedical.com/patents).  
Life moves us is a registered trademark of Globus Medical. Please refer to package insert for  
description, indications, contraindications, warnings, precautions and other important information.

CE 0297

GMTGD187  
12.20 Rev D